

Historic, Archive Document

Do not assume content reflects current scientific knowledge, policies, or practices.

Reserve
aSB998
.N4P58
1957

PLANT NEMATOLOGY
LABORATORY MANUAL

— 0 —

U.S.D.A., NAL

MAR 11 2005

CATALOGING PREP

Prepared by the
Nematology Section

U.S.D.A.

1957

Note: This manual was prepared for use in the
Nematology Section only. It is not for publi-
cation and should not be quoted, cited or copied.

**United States
Department of
Agriculture**



National Agricultural Library

Collecting Soil and Plant Samples

In collecting plant and soil samples for nematode examination, it should be kept in mind that the nematodes feed on roots of growing plants and can be most easily found where their food is abundant (Fig. 1).

A good way to proceed in sampling soil around the roots of growing plants is as follows:

1. Get a shovel, trowel, soil sampling tube or other digging implements.

Also a supply of plastic (or paper) bags, and some labels.

2. Select the plant to sample. The best plants to sample are those which are not yet seriously damaged and still have plenty of living roots. Plants with mostly dead roots will not have many plant parasitic nematodes because there is little food for them.

3. Dig samples to include living roots and some of the adjacent soil. If the plants are small, dig up a whole root system. If the plant is larger, take only part of the root system. A sample of about 500 cc is sufficient.

4. Place the sample in the bag and label it.

5. Protect samples from drying or from excessive heat. Either may kill many of the nematodes.

6. Extract the nematodes from the samples as soon as possible, using the methods described below.

7. If the nematodes are in the plant tissue, such as roots, stems or bulbs, collect plants which show the symptoms and keep from drying until the nematodes are extracted.

Isolations of Nematodes from Soil Samples

Nematodes can be isolated from soil samples by: (1) the Cobb sieving and gravity method and (2) the Baermann funnel.

Sieving and Gravity Method for Soil

1. Put a soil sample of about 300 cc of soil in a bucket (Fig. 2A).
2. Add about 2 liters of water.
3. Mix the soil and water by stirring with a stick, then allow mixture to stand for 30 seconds.
4. Pour the water through a 20-mesh sieve into a second bucket (Fig. 2B).
5. Discard the debris on the 20-mesh sieve.

Note: A larger proportion of the nematodes can be obtained by repeating 2 to 5 once or twice.

6. Discard the material in the first bucket. This consists of sand and heavy soil particles.

7. Pour the water in the second bucket back through a 60-mesh sieve into the first bucket (Fig. 3C).

8. Wash the residue from the 60-mesh sieve into a beaker (Fig. 4E).
9. Pour the water in the first bucket through a 200-mesh sieve.
10. Wash the residue from the 200-mesh sieve into a beaker (Fig. 4E).
11. Pour the material in the beakers into shallow dishes for examination by the microscope. Cysts or adult females of the genus Heterodera will be found in the residue from the 60-mesh sieve. Cysts often float, so both the bottom of the dish and the surface of the water should be examined. Some of the larger nematodes will also be found in the residue from the 60-mesh sieve; the smaller nematodes will be found in the residue from the 200-mesh sieve. Or: Place it in the Baermann funnel as described below.

Isolation of Nematodes from Soil Samples

Nematodes can be isolated from soil samples by (1) the Cobb sieving and gravity method and (2) the Baermann funnel.

Sieving and Gravity Method for Soil

1. Put a soil sample of about 300 cc of soil in a bucket (Fig. 2A).
2. Add about 2 liters of water.
3. Mix the soil and water by stirring with a stick, then allow mixture to stand for 30 seconds.
4. Pour the water through a 20-mesh sieve into a second bucket (Fig. 2B).
5. Discard the debris on the 20-mesh sieve.

Note: A larger proportion of the nematodes can be obtained by repeating 2 to 5 once or twice.

6. Discard the material in the first bucket. This consists of sand and heavy soil particles.
7. Pour the water in the second bucket back through a 60-mesh sieve into the first bucket (Fig. 3C).
8. Wash the residue from the 60-mesh sieve into a beaker (Fig. 4E).
9. Pour the water in the first bucket through a 200-mesh sieve.
10. Wash the residue from the 200-mesh sieve into a beaker (Fig. 4F).
11. Pour the material in the beakers into shallow dishes for examination by the microscope. Cysts or adult females of the genus Heterodera will be found in the residue from the 60-mesh sieve. Cysts often float, so both the bottom of the dish and the surface of the water should be examined. Some of the larger nematodes will also be found in the residue from the 60-mesh sieve; the smaller nematodes will be found in the residue from the 200-mesh sieve. Or: Place it in the Baermann funnel as described below.

Note: Ordinary soil sieves, obtainable from scientific supply companies, can be used for extracting nematodes by the sieving method. The "20-mesh" sieve has 20 openings to the inch and each opening 0.840 mm. wide. The "60-mesh" sieve has openings 0.250 mm. wide and the openings in the "200-mesh" sieve are 0.074 mm. wide.

The Baermann Funnel

A Baermann funnel is made by attaching a short piece of rubber tubing to a funnel and placing a pinchcock on the tubing. The funnel is then placed in an upright position and partly filled with water (Fig. 4F).

1. Place the sample from which the nematodes are to be extracted in a beaker.
2. Place a piece of cloth over the beaker and fasten it with a rubber band.
3. Invert the beaker in the funnel so that the cloth is under the surface of the water.
4. Leave the funnel undisturbed for 3 hours or longer.
5. Take a sample from the funnel by opening the pinchcock for a short time.

The sample should not be more than 10 ml.

6. Place the sample in a shallow dish for examination.

Uses of the Baermann Funnel

The Baermann funnel can be used for extraction of nematodes from (1) soil samples, (2) sieve residues obtained by sieving soil as described above and (3) plant material.

For Soil Samples

1. Mix the soil samples thoroughly.
2. Put about 50 cc in a beaker.
3. Put a cloth cover over the beaker and fasten it with a rubber band.
4. Invert the beaker in the Baermann funnel.

Notes: Ordinary soil sieves, obtainable from scientific supply companies, can be used for extracting nematodes by the sieving method. The "50-mesh" sieve has 50 openings to the inch and each opening 0.850 mm. wide. The "60-mesh" sieve has openings 0.250 mm. wide and the openings in the "100-mesh" sieve are 0.075 mm. wide.

The Baermann Funnel

- A Baermann funnel is made by attaching a short piece of rubber tubing to a funnel and placing a pinchcock on the tubing. The funnel is then placed in an upright position and partly filled with water (Fig. 17).
1. Place the sample from which the nematodes are to be extracted in a beaker.
 2. Place a piece of cloth over the beaker and fasten it with a rubber band.
 3. Invert the beaker in the funnel so that the cloth is under the surface of the water.
 4. Leave the funnel undisturbed for 3 hours or longer.
 5. Take a sample from the funnel by opening the pinchcock for a short time.
- The sample should not be more than 10 ml.

6. Place the sample in a shallow dish for examination.

Uses of the Baermann Funnel

The Baermann funnel can be used for extraction of nematodes from (1) soil samples, (2) sieve residues obtained by sieving soil as described above and (3) plant material.

For Soil Samples

1. Mix the soil samples thoroughly.
2. Put about 50 cc in a beaker.
3. Put a cloth cover over the beaker and fasten it with a rubber band.
4. Invert the beaker in the Baermann funnel.

For Sieve Residues

1. Place a cloth over the beaker containing the sieve residue and fasten it with a rubber band.
2. Invert the beaker in the Baermann funnel.

Plant Material

1. Cut the material into small pieces and place in a beaker.
2. Place a cloth over the beaker and fasten with a rubber band.
3. Invert the beaker in the funnel.

Notes

Extraction of nematodes by the Baermann funnel depends on active movements of the nematodes which work their way through the cloth. Dead nematodes will not be extracted by this method, so it cannot be used with preserved material. In addition to cloth, various kinds of paper have been used in the Baermann funnel. Paper towels work very well as do different kinds of thin porous paper.

Besides the method mentioned, there are numerous ways of supporting the cloth in the Baermann funnel. Wire rings, screen wire or clothes pins have been used.

Nematodes left too long in the Baermann funnel sometimes die from lack of oxygen.

Microscopes and Accessories

The nematologist needs two microscopes as follows:

1. Dissecting Microscope with stand and mirror for substage illumination, equipped with oculars and objectives to give magnifications of about 6 to 10X, 20 to 30X, and 50 to 75X. The lowest magnifications (about 6 to 30X) are used for examination of plant samples for symptoms of nematode attack, and for dissecting plant tissue to isolate the contained nematodes. These magnifications are also used in the various steps of preparing slides. Using the higher magnifications, the skilled worker can make tentative identifications of some kinds of nematodes.
2. Compound microscope equipped with a condensor, three objectives and oculars to give magnifications of about 100X, 450X and 950X. The usual combination is objectives of 16 mm., 4 mm., and about 2 mm. focal lengths, and oculars of about 10X. The lowest magnification is used for locating nematodes on the slide, and for observation of the general size and shape of the body. With magnifications of about 450X, some details of internal and external structure can be seen sufficiently well to permit identification to genus and sometimes to species. The highest magnification, obtained with an oil immersion objective, is used for study of minute details and is necessary for identification of most species.

A light which can be used for either incident or transmitted illumination is needed for the dissecting microscope, and a good microscope lamp is needed for the compound microscope.

An eyepiece micrometer or filar micrometer is a useful accessory for the compound microscope.

A camera lucida is needed if drawings of nematodes are to be made.

The nematologist needs two microscopes as follows:

1. Dissecting Microscope with stand and mirror for oblique illumination, equipped with oculars and objectives to give magnifications of about 6 to 10X, 20 to 30X, and 50 to 75X. The lowest magnifications (about 6 to 10X) are used for examination of plant samples for symptoms of nematode attack, and for dissecting plant tissue to isolate the contained nematodes. These magnifications are also used in the various steps of preparing slides. Using the higher magnifications, the skilled worker can make tentative identifications of some kinds of nematodes.
 2. Compound microscope equipped with a condenser, three objectives and oculars to give magnifications of about 100X, 450X and 950X. The usual combination is objectives of 16 mm., 4 mm., and about 2 mm. focal lengths, and oculars of about 10X. The lowest magnification is used for locating nematodes on the slide, and for observation of the general size and shape of the body. With magnifications of about 450X, some details of internal and external structure can be seen sufficiently well to permit identification to genus and sometimes to species. The highest magnification, obtained with an oil immersion objective, is used for study of minute details and is necessary for identification of most species.
- A light which can be used for either incident or transmitted illumination is needed for the dissecting microscope, and a good microscope lamp is needed for the compound microscope.
- An eyepiece micrometer or filar micrometer is a useful accessory for the compound microscope.
- A camera lucida is needed if drawings of nematodes are to be made.

Finding Nematodes in Plant Tissue

Nematodes in plant tissue may be (1) active nematodes or (2) inactive. Different methods are used in searching for each.

Active Nematodes

Active nematodes are able to move and are always slender. This kind of nematode can be found in plant tissue as follows:

1. Place a small piece of the plant tissue in a watch glass and cover with water (Fig. 5A).
2. Place the watch glass under the dissecting microscope, using a magnification of 10x to 30x.
3. With dissecting needles, tear the tissue apart (Fig. 6B).
4. This frees the nematodes from the plant tissue and they can then be found in the water. Look for them in the bottom of the dish.
5. If no nematodes are found, set the dish aside for an hour or more. Then look again.

A very efficient method of finding active nematodes in roots is the "Incubation Technique" described by Young in 1954 (Fig. 7).

1. Wash the roots thoroughly.
2. Place the wet roots in a fruit jar with a cover.
3. Set the jar aside for several days.
4. Spray a small amount of water on the roots and then pour it off into a beaker.
5. Examine this water under the dissecting microscope.

The Waring Blendor

When nematodes are numerous in plant tissue, they can be found as follows:

1. Place about 2 grams of washed plant material in a Waring Blendor and add about 100 ml. of water (Fig. 8A).
2. Run the Blendor for 20 seconds.

3. Pour the material through a 60-mesh sieve over a 200-mesh sieve (Fig. 8B).
4. Wash the material on the sieves with a gentle stream of tap water.
5. Discard the residue on the 60-mesh sieve and look for nematodes in the residue on the 200-mesh sieve.
6. Or: Place the residue from the 200-mesh sieve in the Baermann funnel.

Inactive Nematodes

Inactive nematodes are the species with "pear-shaped" or otherwise enlarged bodies, such as the root-knot nematodes. These will be partially or completely embedded in the plant tissue.

1. Place short pieces of the plant tissue in a watch glass under the dissecting microscope.
2. Using dissecting needles, pull the tissue apart to expose the nematodes.
3. Dissect the nematodes out of the plant tissue.

Note: Root-knot nematodes will be found mostly in roots with distinct knots, but are occasionally found in roots without distinct knots. The nematodes are in the knots, and their exact location is often marked by the presence of an egg mass. This is an irregular brown body about one millimeter in diameter on the surface of the root. It can easily be detached, revealing the nematode underneath.

Root-knot nematodes are more easily dissected out of roots which have been in 5% formalin for at least 24 hours than from fresh roots.

Mounting Nematodes on Temporary Microscope Slides

Microscope slides of nematodes which will be usable for several weeks or months can be made as follows:

1. Place a drop of 5% formaldehyde in the center of a clean microscope slide (Fig. 10A).
2. Locate a nematode in the Syracuse Watch Glass with the dissecting microscope, using a magnification of about 15 to 30X.
3. Gently lift the nematode with the point of the bamboo splinter upward to the water surface. Keep it in focus under the microscope (Fig. 9).
4. When the nematode is just under the surface of the water, lift rapidly to bring it through the surface film.
5. Place the point of the bamboo splinter in the drop of 5% formaldehyde on the slide and withdraw it carefully. This will usually leave the nematodes in the drop (Fig. 10A).
6. Place the slide under the dissecting microscope and make sure that all the nematodes are at the bottom of the drop.
7. Using forceps, carefully place a cover glass over the drop (Fig. 10B).
8. Still working under the microscope, use a small piece of filter paper, blotting paper or paper towel to absorb the excess liquid. Be careful that none of the nematodes are lost (Fig. 11A).
9. Seal the slide with either: (a) A mixture of one part of paraffin and one part of vaseline. This is kept hot and applied with a small brush. (b) ZUT, a special slide-sealing compound obtainable from Bennett's, 65 West First South Street, Salt Lake City 10, Utah. (c) Other materials including fingernail polish (Fig. 11B).
10. Label the slide (Fig. 11C).

Microscopic slides of nematodes which will be ready for several weeks or

months can be made as follows:

1. Place a drop of 5% formaldehyde in the center of a clean microscope

slide (Fig. 10A).

2. Locate a nematode in the Syracuse Watch Glass with the dissecting

microscope, using a magnification of about 15 to 30X.

3. Gently lift the nematode with the point of the bamboo splinter upward

to the water surface. Keep it in focus under the microscope (Fig. 9).

4. When the nematode is just under the surface of the water, lift rapidly to

bring it through the surface film.

5. Place the point of the bamboo splinter in the drop of 5% formaldehyde on

the slide and withdraw it carefully. This will usually leave the nematode in

the drop (Fig. 10A).

6. Place the slide under the dissecting microscope and make sure that all

the nematodes are at the bottom of the drop.

7. Using forceps, carefully place a cover glass over the drop (Fig. 10B).

8. Still working under the microscope, use a small piece of filter paper,

blotting paper or paper towel to absorb the excess liquid. Be careful that none

of the nematodes are lost (Fig. 11A).

9. Seal the slide with either: (a) A mixture of one part of paraffin and

one part of vaseline. This is kept hot and applied with a small brush. (b) BUT,

a special slide-sealing compound obtainable from Bennett's, 65 West First South

Street, Salt Lake City 10, Utah. (c) Other materials including fingernail

polish (Fig. 11B).

10. Label the slide (Fig. 11C).

Notes: A bamboo splinter suitable for picking up nematodes is made by splitting a piece of bamboo about 20 cm long to obtain a piece about 2 mm wide. One end of this ^{is}/sharpened with a knife to make a long and very fine point. The point is best made by trimming with a razor blade under the dissecting microscope. The finer the point, the easier it is to use (Fig. 9).

Cover glasses should be the thinnest obtainable, usually the "No. 1" or "No. 0" thickness. If cover glasses are too thick, the nematodes cannot be examined with the oil immersion lens. Diameter of cover glasses should be 12 to 18 mm.

The slide sealing mixture is painted on, working under the dissecting microscope to be sure that the edges of the cover slip are sealed.

Other Methods of Making Microscope Slides of Nematodes

1. Temporary slides for immediate use under the lower powers of the compound microscope can be made by simply placing the nematodes in a drop of water and then placing a cover glass. If the nematodes are moving, heat the slide very gently with a match. Try to heat just enough to stop movement.

2. Slides can also be made by mounting nematodes which have been preserved in 5% formaldehyde in lacto-phenol solution. (Melted phenol, 3 parts; lactic acid, 1 part; glycerine, 2 parts; water, 1 part.) Such slides deteriorate in a few months so that certain internal structures of the nematode are difficult to see. Cuticle details remain good. This is an excellent mounting medium for cysts of Heterodera species and for perineal patterns of Meloidogyne species.

3. Permanent slides are made by mounting nematodes in glycerine. See below under "Fixing and Preserving Specimens".

Fixing and Preserving Specimens

The standard material for preserving nematodes or plants and soil containing nematodes is 5% formaldehyde. A good procedure is as follows:

1. Prepare 10% formaldehyde by proper dilution with water of commercial formaldehyde solutions containing about 36% to 40% formaldehyde.

2. Place the material to be fixed in water at room temperature or below and add an equal volume of 10% formaldehyde solution heated to just below the boiling point. This relaxes and kills the nematodes without danger of ruining them by overheating.

3. Store in tightly closed containers until ready for examination. Or, transfer to pure glycerine as described below.

Note: Nematodes, or plant material containing nematodes will remain in good condition for years in 5% formaldehyde. Most other materials which have been tried have given inferior results in the long run. Formalin can also be used cold and this is sometimes an advantage if cuticular structures are to be examined. If sufficient material is at hand, it may be desirable to divide it and use both hot and cold formaldehyde.

4. Prepare solution of 1.5% glycerine in 10% ethyl alcohol.

5. Place the nematodes in this solution in Bureau of Plant Industry Watch Glasses or other small shallow containers. (Note B.P.I. Watch Glasses are Syracuse watch glasses 27 mm. in diameter.)

6. Place the watch glass in a tightly closed container of about 200 ml. capacity with a 16 mm. by 57 mm. screw cap vial filled with desiccated calcium carbonate. The cap of the vial should have a hole about 2 mm. in diameter. The purpose of this procedure is to permit very gradual evaporation of the water and alcohol, leaving the nematodes in pure glycerine after 3 to 6 weeks. If evaporation is too rapid, the nematodes will often collapse.

material for reservoirs of plants and soil containing
 formaldehyde. A good procedure is as follows:

1. From 10% formaldehyde by proper dilution with water of commercial formaldehyde solutions containing about 36% to 40% formaldehyde.
2. Place the material to be fixed in water at room temperature or below and add an equal volume of 10% formaldehyde solution heated to just below boiling point. This relaxes and kills the nematodes without danger of retaining by overheating.

Store in tightly closed containers until ready for examination.

Material to pure glycerine as described below.

Materials: Nematodes, or plant material containing nematodes will remain in good condition for years in 5% formaldehyde. Most other materials which have been tried give inferior results in the long run. Formalin can also be used cold and will give an advantage if cylindrical structures are to be examined. If material is at hand, it may be desirable to divide it and use both hot and cold formaldehyde.

1. Solution of 1.5% glycerine in 10% ethyl alcohol.

2. Place the nematodes in this solution in Bureau of Plant Industry Watch

glass or other small shallow containers. (Note B.P.I. Watch Glasses are

available in glasses 27 mm. in diameter.)

3. Place the watch glass in a tightly closed container of about 100 ml.

4. Fill with a 10 mm. by 25 mm. screw cap vial filled with deaerated calcium

hydroxide. The cap of the vial should have a hole about 2 mm. in diameter. The

purpose of this procedure is to prevent very gradual evaporation of the water and

to keep the nematodes in pure glycerine for 3 to 6 weeks. If evaporation

7. Transfer the watch glasses to a large desiccator containing calcium chloride until it is certain that all the water and alcohol have been evaporated. The nematodes can then be stored indefinitely in desiccated glycerine in tightly closed containers.

Making Permanent Slides

1. Keep a shallow dish of glycerin for use as the mounting medium and another dish containing short pieces of glass rod in the desiccator. The glass rods should have about the same diameter as nematodes, that is, should be of varying diameters in the range of 10 to 150 microns. These rods are obtained by cutting glass wool into lengths of several millimeters, or by drawing out threads from glass rods heated by a Bunsen burner.

2. Place a small drop of the mounting medium on a slide and transfer nematodes to it from the watch glasses.

3. Select three pieces of glass rod just a little smaller in diameter than the nematodes and place at the edges of the drops evenly spaced. These are to support the cover glass so that the nematodes are not flattened or distorted.

4. Lower a cover glass gently over the drop of mounting fluid.

5. Blot up excess glycerine from around the edges of the cover glass.

6. Seal the slide with ZUT.

Note: Nematodes can be mounted on ordinary microscope slides, but there is a considerable advantage in the use of the metal slides invented by Cobb (1917). These consist of a piece of aluminum bent to hold a square cover glass and 2 pieces of cardboard. Nematodes are mounted between the square cover glass and a round cover glass and so can be examined from either side. The cardboard pieces are used as labels. Metal slides can be stacked and so stored in boxes without the slots which are needed for glass slides.

7. Transfer the watch glasses to a large desiccator containing calcium chloride until it is certain that all the water and alcohol have been evaporated. The desiccator can then be stored indefinitely in desiccated glycerine in tightly sealed containers.

8. Keep a shallow dish of glycerine for use as the mounting medium and another dish containing short pieces of glass rod in the desiccator. The glass rods should have about the same diameter as nematodes, that is, should be of uniform diameter in the range of 10 to 150 microns. These rods are obtained by cutting glass wool into lengths of several millimeters, or by drawing out threads from glass rods heated by a Bunsen burner.

9. Place a small drop of the mounting medium on a slide and transfer

nematodes to it from the watch glasses.

10. Select three pieces of glass rod just a little smaller in diameter than the nematodes and place at the edges of the drops evenly spaced. These are to support the cover glass so that the nematodes are not flattened or distorted.

11. Lower a cover glass gently over the drop of mounting fluid.

12. Blot up excess glycerine from around the edges of the cover glass.

13. Seal the slide with EUT.

Note: Nematodes can be mounted on ordinary microscope slides, but there is

an undeniable advantage in the use of the metal slides invented by Cobb (1917).

These consist of a piece of aluminum bent to hold a square cover glass and 2 pieces

of cover glass. Nematodes are mounted between the square cover glass and a round

cover glass and so can be examined from either side. The cardboard pieces are

used as spacers. Metal slides can be stacked and so stored in boxes without the

at the ends of the slides for glass slides.

Nematode Anatomy

All nematodes at some stage of their life are elongate worms, the body shape ranging from nearly cylindrical with rounded ends to fusiform. Most nematodes are this shape all through their lives, but adult females of certain important plant parasites have enlarged bodies which are pear-shaped, lemon-shaped, kidney-shaped, or otherwise enlarged (Fig. 12). The only appendages nematodes have are short setae usually near the anterior end. Soil and plant parasitic nematodes are rather small, average length being about a millimeter, with a minimum length of about 0.25 mm. and a maximum of about 10 mm. The great majority of the species are less than 2 mm. long.

In the typical nematodes (Fig. 14) the mouth is terminal at the anterior end. There is a mouth cavity ^{is} which/reduced to a stylet in plant parasitic forms. Attached to the mouth cavity or stylet is the oesophagus which leads to the intestine. The intestine leads to the rectum and ends at the anus, which is situated on the ventral side. The body from the anus to the posterior terminus is the tail.

In the female, the vulva is also situated on the ventral side and is a transverse slit. This leads to the ovaries which are tubes lying in the body cavity. There may be either one or two ovaries. At the posterior end of the male body are a pair of structures known as spicules. These normally lie in the body, but can be pushed out through the anus. They are used in copulation. Also leading to the anus of the male is a single tube in which sperms are formed called the testis. Rarely, there are two testes.

Nematodes also have an excretory system opening into a ventral excretory pore situated on the anterior part of the body, in most species. The excretory canal, leading to this pore is the only part of this system which can be seen.

Nematode Anatomy

All nematodes at some stage of their life are elongate worms, the body shape varying from nearly cylindrical with rounded ends to fusiform. Most nematodes are of a shape all through their lives, but adult females of certain important plant parasites have enlarged bodies which are pear-shaped, lemon-shaped, kidney-shaped, or otherwise enlarged (Fig. 12). The only appendages nematodes have are short bristles usually near the anterior end. Soil and plant parasitic nematodes are small, average length being about a millimeter, with a minimum length of about 0.25 mm. and a maximum of about 10 mm. The great majority of the species are less than 2 mm. long.

In the typical nematodes (Fig. 11) the mouth is terminal at the anterior end. There is a mouth cavity which is reduced to a stylet in plant parasitic forms. Attached to the mouth cavity or stylet is the oesophagus which leads to the intestine. The intestine leads to the rectum and ends at the anus, which is situated on the ventral side. The body from the anus to the posterior terminus is the tail.

In the female, the vulva is also situated on the ventral side and is a transverse slit. This leads to the ovaries which are tubes lying in the body cavity. There may be either one or two ovaries. At the posterior end of the male there are a pair of structures known as spermathecae. These normally lie in the body, but can be pushed out through the anus. They are used in copulation. Also leading to the anus of the male is a single tube in which sperms are formed. In the testis, however, there are two testes.

Nematodes also have an excretory system opening into a ventral excretory pore situated in the anterior part of the body, in most species. The excretory pore, leading to this pore is the only part of this system which can be seen.

There is also a nervous system, the only visible part of which is the nerve ring around the oesophagus, and this is usually difficult to see.

The body is lined with muscle fibers arranged in four sectors separated by two lateral chords, a dorsal and a ventral chord. These are usually difficult to see.

Sensory organs include papillae surrounding the mouth cavity, near the excretory pore and near the posterior end. At the anterior end are also two chemical sense organs called amphids. In the Phasmidia these are terminal and very small so are seldom seen. In certain of the Aphasmidia, the external amphids are conspicuous circles, spirals and similar forms.

The phasmids are generally located in the tail. These are two glands connecting with lateral pores located in the middle of the lateral fields. Usually only the pores can be seen and often these are hard to locate. The canals leading to the surface can sometimes be seen in ventral or dorsal views.

The body of the nematode is covered with a transparent, colorless or slightly yellowish cuticle. This cuticle may be marked in various ways. Many nematodes have annules running at right angles to the body axis. These may be less than one micron wide or as much as 5 microns wide. Often they are interrupted laterally by the "lateral fields" which are more or less elevated ridges running lengthwise of the nematode. In many Tylenchida, the lateral fields are divided into three or more strips by lines.

In one genus of the Tylenchidae (Criconema) the coarse annules are bordered by well developed rows of spines or scales.

There is also a nervous system, the only visible part of which is the nerve

running around the esophagus, and this is usually difficult to see.

The body is lined with muscle fibers arranged in four sections separated by

two lateral chords, a dorsal and a ventral chord. These are usually difficult

to see.

There are also organs (mouth papillae) surrounding the mouth cavity, near the

esophageal pore and near the posterior end. At the anterior end are also two

chemical sense organs called ampullae. In the Trematoda these are ventral and

very small; in the Agnatha, the external

ampullae are conspicuous circles, spirals and similar forms.

The pharynx is generally located in the tail. There are two glands

connecting with lateral pores located in the middle of the lateral fields.

Only the pores can be seen and often these are hard to locate. The

canals leading to the surface can sometimes be seen in ventral or dorsal views.

The body of the nematode is covered with a transparent, colorless or

slightly yellowish cuticle. This cuticle may be marked in various ways. Many

nematodes have annules running at right angles to the body axis. These may be

less than one micron wide or as much as 5 microns wide. Often they are interrupted

laterally by the "lateral fields" which are more or less elevated ridges running

longitudinally as the nematode. In many Trematoda, the lateral fields are divided

into dorsal and ventral parts.

The body of the nematode is usually cylindrical, but may be flattened or

well developed. The head is usually at one end.

There is a small mouth at the anterior end.

The body is usually covered with a thin cuticle.

The body is usually covered with a thin cuticle.

The body is usually covered with a thin cuticle.

Nematode Systematics

Nematodes belong to the Phylum Nematoda according to some authors, and to the Phylum Aschelmintha according to others. It is generally agreed that they should be subdivided into 2 classes as suggested by Chitwood in 1933. These classes are (1) Phasmidia, and (2) Aphasmidia; that is, nematodes with phasmids and nematodes without phasmids.

The classes are divided into orders with names derived from that of the type genus and ending in "-ida" (example Rhabditida). Names of suborders end in "-ina" (Rhabditina); of superfamilies in "-oidea" (Rhabditoidea); of families in "-idae" (Rhabditidae); and of subfamilies in "-inae" (Rhabditinae). The families are composed of genera and the genera of species. A few subspecies have been named.

About 5000 nematode species which are parasites of vertebrates are known. About an equal number of species of soil, marine, freshwater, insect parasites and plant parasites are known. This may be only a tenth or less of the existing species. The nematode fauna of the soil and plants of the greater part of the world is unknown.

In this manual we are concerned principally with the plant parasitic nematodes. Since the plant nematologist often encounters free-living nematodes in his work, we will also discuss the most common genera of these.

The plant parasitic forms now known belong mostly to a single order, the Tylenchida, of the class Phasmidia. A few species of 3 genera of the superfamily Dorylaimoidea of the Aphasmidia are also plant parasites. Soil nematodes are found in both Phasmidia and Aphasmidia.

The best book on taxonomy of the soil and plant parasitic nematodes is "Soil and Freshwater Nematodes", by T. Goodey. This was published in 1951 by Methuen & Co., Ltd., London and by John Wiley & Sons, Inc., New York. It contains description of the genera with which we are concerned and illustrations

of a representative species of each genus. The next best book is "A Manual of Agricultural Helminthology", by I. N. Filipjev and J. H. Schuurmans-Stekhoven, published by E. J. Brill, Leiden, Holland in 1941. This is obtainable from E. G. Stechert & Co., 31 East 10th Street, New York City.

With the keys in this manual and "Soil and Freshwater Nematodes", identification of the common soil and plant parasitic nematodes to genus is possible.

Nematode taxonomy is constantly growing and changing. New species of nematodes are constantly being found, and new genera are fairly common. In addition a certain amount of changing of species from one genus to another is constantly being done. Consequently, identification of nematodes to species is difficult and often nearly impossible unless the genus has been recently revised.

In using the keys, constant reference to the illustrations of this manual and to those of the reference books should be made. There is no substitute for illustrations.

Identification of Nematodes

The attached keys are used for nematode identification. The "Key to the most Common Nematodes of Agricultural Soils and Plants" will facilitate identification to family or subfamily.

As the first step a decision is made as to whether the nematode belongs to the Phasmidia or Aphasmidia. This is best determined by the form of the oesophagus. In the Phasmidia we find the oesophagi like those shown in figures 15 and 16.

Note that the "Rhabditoid" oesophagus has a distinct bulb at the base and that the bulb contains a valve of distinctive shape. This valve is easy to see (Fig. 15A).

The "Diplogasteroid" oesophagus is of the same general shape but the valve is lacking (Fig. 15B).

The other 3 types of oesophagi, "Tylenchoid", "Neotylenchoid" and "Aphelenchoid" are characterized by a stylet which is always present and easy to see. Stylets may be variations of any of the forms shown in figure 17 A-E. Attached to the stylet is a thin oesophageal tube. In all except the Neotylenchoid oesophagus (Fig. 16C), this tube leads to an oval "Valve" in a "median oesophageal bulb", then continues on to the intestine.

Males of some species of the Phasmidia may have a distinct bursa as shown in figure 14 and in figure 19D. Mature females of some genera have enlarged bodies as shown in figure 12.

In the Aphasmidia, the most common types of oesophagi are the "Plectoid", "Cylindrical" and "Dorylaimoid" as shown in figure 17A-C.

Note that the "Plectoid" oesophagus has a basal bulb as in the Rhabditoid type, but that the valve is elongated (Fig. 17C).

Identification of Nematodes

The attached keys are used for nematode identification. The "Key to the Common Nematodes of Agricultural Soils and Plants" will facilitate identification for the family or subfamily.

As the first step a decision is made as to whether the nematode belongs to the Phasmidae or Aphasmidae. This is best determined by the form of the oesophagus.

In the Phasmidae we find the oesophagi like those shown in figures 15 and 16.

Note that the "Phasmid" oesophagus has a distinct bulb at the base and that the bulb contains a valve of distinctive shape. This valve is easy to see

(fig. 15A).

The "Diplogasteroid" oesophagus is of the same general shape but the valve

is lacking (fig. 15B).

The other 3 types of oesophagi, "Neostylenoid", "Stylenoid", and "Aphasmidoid"

are characterized by a stylet which is always present and easy to see. Stylets may

be located at any of the forms shown in figure 17 A-E. Attached to the stylet

is a thin oesophageal tube. In all except the Neostylenoid oesophagus (fig. 16C),

this tube leads to an oval "Valve" in a "median oesophageal bulb", then continues

on to the intestine.

Most of some species of the Phasmidae may have a distinct bursa as shown

in figure 18 and in figure 19D. Mature females of some genera have enlarged bodies

as shown in figure 18.

In the Aphasmidae, the most common types of oesophagi are the "Plicatoid",

"Cylindrical" and "Dorylaimoid" as shown in figure 17A-C.

Note that the "Plicatoid" oesophagus has a basal bulb as in the Phasmidoid

type, but that the valve is elongated (fig. 17C).

The cylindrical type of oesophagus is easy to recognize. It is seen most often in the genus Mononchus, where it is associated with a large mouth cavity containing a tooth as shown in figure 17A.

The "Dorylaimoid" type of oesophagus consists of a narrow anterior portion and a broader posterior portion as shown in figure 17B. This is often, though not always associated with stylets of the types shown in figure 18F-H.

If the nematode belongs to the Phasmidia and has a stylet it will be placed in the Tylenchida, as shown by Item 2 of the Key. If the stylet is absent, it belongs to the Rhabditida.

Nematodes belonging to the Tylenchida are separated into Tylenchoidea and Aphelenchoidea by the characters named in Item 3 of the Key. If the stylet has well developed knobs as shown in figure 18A-D, the nematode belongs to the Tylenchoidea. The oesophageal tube in this group usually has an abrupt bend just posterior to the stylet. At this point can be seen a short branch which is the opening of the dorsal oesophageal gland. If this is present and the stylet knobs are absent, the nematode should have a long tail as shown in figure 19E-F.

On the other hand, if the nematode has a large prominent oesophageal bulb and a stylet with very small knobs or no knobs, it belongs in the Aphelenchoidea (Fig. 16B).

Phasmidia without stylets are placed in the Diplogasteridae if there is no valve in the basal oesophageal bulb (Fig. 15B).

If there is a valve in the dorsal oesophageal bulb, and the stoma is cylindrical as shown in figure 15A, the nematode is placed in the Rhabditoidea. If the stoma is not cylindrical, the nematode is placed in the Cephalobidae.

Nematodes belonging to the Aphasmidia are placed in the plectinae if the oesophagus is as shown in figure 18C. Most nematodes of the group have a spinneret on the tail as shown in figure 19H.

The cylindrical type of oesophagus is easy to recognize. It is seen most often in the *Rhabditidae*, where it is associated with a large mouth cavity (see figure 17A).

The "bottle-neck" type of oesophagus consists of a narrow anterior portion and a broader posterior portion as shown in figure 17B. This is often, though not always associated with stylets of the types shown in figure 18B-H.

If the nematode belongs to the *Rhabditidae* and has a stylet it will be placed in the *Tylenchidae*, as shown by Item 2 of the Key. If the stylet is absent, it belongs to the *Rhabditidae*.

Nematodes belonging to the *Tylenchidae* are separated into *Tylenchidae* and

des by the characters named in Item 3 of the Key. If the stylet

has well developed knobs as shown in figure 18A-D, the nematode belongs to the

Tylenchidae. The oesophageal tube in this group usually has an abrupt bend

just posterior to the stylet. At this point can be seen a short branch which

is the opening of the dorsal oesophageal gland. If this is present and the

stylet knobs are absent, the nematode should have a long tail as shown in figure

On the other hand, if the nematode has a large prominent oesophageal bulb

and a stylet with very small knobs or no knobs, it belongs in the *Apeloneuridae*

(fig. 18E).

Rhabditidae without stylets are placed in the *Allogasteridae* if there is no

in the basal oesophageal bulb (fig. 18F).

If there is a valve in the dorsal oesophageal bulb, and the stoma is

cylindrical as shown in figure 18A, the nematode is placed in the *Rhabditidae*.

If the stoma is not cylindrical, the nematode is placed in the *Tylenchidae*.

Nematodes belonging to the *Apeloneuridae* are placed in the *Apeloneuridae* if the

oesophagus is as shown in figure 18C. Most nematodes of the group have a

on the tail as shown in figure 18H.

If the oesophagus is cylindrical, as shown in figure 17A, the nematode is placed in the Mononchidae.

If the oesophagus is dorylaimoid, as shown in figure 17B, the nematode is placed in the Dorylaimoidea.

To identify a nematode:

1. Study the oesophagus. Look at four or five specimens if possible.
2. Compare the oesophagus with figures 15, 16 and 17, and decide which type it is.
3. Use the key to identify the nematode to superfamily or family.
4. Find the section relating to this group in "Soil and Freshwater Nematodes".
5. Compare your nematode with the illustrations in the book. This will usually result in an identification to genus.

...the nematode is ... in figure 17A, the nematode is ...

...the nematode is ... in figure 17B, the nematode is ...

...the nematode is ... in figure 17C, the nematode is ...

...the nematode is ... in figure 17D, the nematode is ...

...the nematode is ... in figure 17E, the nematode is ...

...the nematode is ... in figure 17F, the nematode is ...

...the nematode is ... in figure 17G, the nematode is ...

...the nematode is ... in figure 17H, the nematode is ...

...the nematode is ... in figure 17I, the nematode is ...

...the nematode is ... in figure 17J, the nematode is ...

...the nematode is ... in figure 17K, the nematode is ...

...the nematode is ... in figure 17L, the nematode is ...

...the nematode is ... in figure 17M, the nematode is ...

...the nematode is ... in figure 17N, the nematode is ...

...the nematode is ... in figure 17O, the nematode is ...

...the nematode is ... in figure 17P, the nematode is ...

...the nematode is ... in figure 17Q, the nematode is ...

...the nematode is ... in figure 17R, the nematode is ...

...the nematode is ... in figure 17S, the nematode is ...

...the nematode is ... in figure 17T, the nematode is ...

...the nematode is ... in figure 17U, the nematode is ...

...the nematode is ... in figure 17V, the nematode is ...

...the nematode is ... in figure 17W, the nematode is ...

...the nematode is ... in figure 17X, the nematode is ...

...the nematode is ... in figure 17Y, the nematode is ...

...the nematode is ... in figure 17Z, the nematode is ...

Key to the Most Common Nematodes of Agricultural Soils & Plants

1. Oesophagus rhabditoid (with valve in basal bulb as in Fig. 15A),
diplogasteroid (without valve in basal bulb as in Fig. 15B),
tylenchoid (Fig. 16A & B), or aphelenchoid (Fig. 16C). If
tylenchoid or aphelenchoid, always with stylet. Males of some
genera have a distinct bursa (Fig. 19D). Phasmids always present
but most often difficult to locate. Mature females of some genera
have a much enlarged body Phasmidia 2.
- Oesophagus plectoid (Fig. 17C), cylindrical (Fig. 17A) or dorylaimoid
(Fig. 17B). Often with setae. Males without bursa. Always with
elongate, vermiform body Aphasmidia 6.
2. Stylet always present. Shapes and proportions of stylets vary
greatly with genera and species, but are usually recognizable
as variations of the forms shown in Fig. 17A-E. Attached to
the stylet is a thin oesophageal tube which may be straight or
coiled Tylenchida 3.
- Stylet absent. Anterior portion of oesophagus muscular (striated)
(Fig. 15A-B) Rhabditida 4.
3. Stylet mostly with well-developed knobs. If stylet knobs are absent,
tail is long and thin (Fig. 19E & F). Dorsal oesophageal gland
orifice near base of stylet, or at most, not more than one stylet
length posterior to stylet knobs. The oesophageal tube often has
an abrupt bend at this point. (Fig. 16A & B) Tylenchoidea
- Stylet knobs small or absent. Median oesophageal bulb occupying nearly
full width of body (Fig. 16B). Tail rounded or conical (Fig. 19A-B).
Dorsal oesophageal gland orifice in median bulb and difficult to locate.
Oesophageal tube without abrupt bends Aphelenchoidea

1. *Oesophagus* rhachis (with valve in basal bulb as in fig. 12A).

apophagoid (without valve in basal bulb as in fig. 12B).

typhlocyba (fig. 10A-B), or *apophagoid* (fig. 10C). If

typhlocyba or *apophagoid*, always with stylet. Males of some

genera have a distinct bursa (fig. 10D). *Phanidia* always present

but most often difficult to locate. Mature females of some genera

have a much larger body *Phanidia* 2.

Oesophagus pharynx (fig. 11C), cylindrical (fig. 11A) or dorsoventral

(fig. 11B). Often with setae. Males without bursa. Always with

Stylet always present. Shapes and proportions of stylets vary

greatly with genera and species, but are usually recognizable

as variations of the forms shown in fig. 11A-E. Attached to

the stylet is a thin oesophageal tube which may be straight or

curved *Phanidia* 2.

Stylet absent. Anterior portion of oesophagus muscular (stretched)

. *Phanidia* 1.

Stylet mostly with well-developed knobs. If stylet knobs are absent,

tail is long and thin (fig. 12C-E). Dorsal oesophageal gland

orifice near base of stylet, or at most, not more than one stylet

length posterior to stylet knobs. The oesophageal tube often has

an abrupt bend at this point. *Phanidia* 1.

Stylet knobs small or absent. Median oesophageal bulb occupying nearly

full width of body (fig. 12B). Tail rounded or conical (fig. 12A-B).

Dorsal oesophageal gland orifice in median bulb and difficult to locate.

Oesophageal tube without abrupt bends *Phanidia* 1.

4. Oesophagus rhabditoid, with valve in basal bulb (Fig. 15A). 5.

Oesophagus diplogasteroid, without valve in basal bulb (Fig. 15B).

. Diplogasteridae

5. Stoma cylindrical, usually much longer than wide (Fig. 15A).

. Rhabditidae

Stoma not cylindrical, or if nearly so, about as long as wide. Lips

of some genera have distinct projections ranging in shape from

rounded to elaborately ornamented. Cephalobidae

6. Oesophagus plectoid (with basal bulb, Fig. 17C). Tail tip with

a small projection (spinneret) Fig. 19H. Plectinae

Oesophagus cylindroid (Fig. 17A) or dorylaimoid (Fig. 17B). 7.

7. Oesophagus cylindroid. Mouth cavity large, subglobular, usually

with one or more large teeth (Fig. 17A). Mononchidae

Oesophagus dorylaimoid (Fig. 17B). Stylet often as shown in Fig. 13C.

Some with much longer stylet, Fig. 18F-H, others with a tooth. . . Dorylaimoidea

Note: Known plant parasites are either Tylenchida, or species of the

genera Xiphinema, Longidorus and Trichodorus of the Dorylaimoidea.

IDENTIFICATION OF PLANT PARASITIC NEMATODES *

Nematodes which are plant parasites always have a stylet which can be recognized as some variation of the stylet shapes shown in Figure 18. Nematodes which do not have stylets, or which have stylets which are not obvious variations of the shapes shown in Figure 18 are not plant parasites. Having determined that the stylet resembles those shown in Figure 18, the nematode can be identified to genus by use of the key.

In using the key, constant reference should be made to the illustrations. When a tentative identification has been made, the nematode should be compared with the illustrations in "Soil and Freshwater Nematodes", by T. Goodey (published in 1951 by Methuen & Co., Ltd., London, and John Wiley & Sons, Inc., New York). This book illustrates only one species of each genus but other species of the genus will be very similar in most respects.

If a stylet is present, it is easily visible on good specimens at magnification of about 400 times. It is always good practice to examine several specimens; characters may be difficult to see on one specimen and easy to see on another. It should also be noted that the key can be used only for the identification of mature females. This is the stage most frequently encountered; for some of the species males are very rare.

* As used here the term "plant parasites" includes parasites of fungi. It should be kept in mind that many kinds of nematodes probably feed only on fungi, others may feed either on fungi or higher plants. Probably most of the Neotylenchidae feed only on fungi, as do most species of Aphelenchus. Certain species of Aphelenchoides and Ditylenchus are known to feed either on fungi or higher plants, and the same may be true for species of other genera of the Tylenchida as well. Thus, nematodes with stylets as shown in Figure 18 may or may not be parasites of higher plants.

Key to the Mature Females of the Common Plant Parasitic Nematodes

1. With median oesophageal bulb (Figs. 16A & C) 5
 Without median oesophageal bulb as in figure 16B; Or without median
 oesophageal bulb and with stylet as in Fig. 18F, G or H. 2
2. Stylet short (Fig. 16B) Neotylenchidae
 Stylet long (Fig. 18F-H) 3
3. Stylet straight without enlargement at base (Fig. 18F). Longidorus
 Stylet with enlargement at base 4
4. Stylet straight with oval enlargements at base (Fig. 18G). Xiphinema
 Stylet curved (Fig. 18H) Trichodorus
5. Body of mature female pear-shaped, lemon-shaped or enlarged and saccate.
 (Fig. 12A-G.) Found in roots of plants, either embedded or attached
 by neck, some as cysts in soil 6
 Body of mature female much longer than wide (Fig. 12 H-L) 10
6. Body of mature female pear-or lemon-shaped (Fig. 12A-C) 7
 Body of mature female saccate (Fig. 12D-G) 8
7. Body of mature female pear-shaped, white, (Fig. 12A). Found completely
 embedded in roots, nearly always in distinct knots . . . Meloidogyne
 Body of mature female pear-shaped or lemon-shaped, (Fig. 12B-C) white,
 yellowish, or brown according to age, attached to root by neck only,
 or found as a brown cyst in soil Heterodera
8. Mature female embedded in plant root, often in knot, body shape ovoid to
 spheroid with elongated "tail", vulva nearly terminal.
 (Fig. 12F-G) Nacobbus
 Mature females attached to root by neck, body more or less kidney-shaped.
 (Fig. 12D-E) 9
9. Vulva at 90% of body length, (Fig. 12E) parasites of citrus and olives . .
 Tylenchulus
 Vulva at 72% of body length, parasites of numerous plants, mostly annuals.
 (Fig. 12D) Rotylenchulus

10. Tail long and thin, more than 6 times anal body diameter (Fig. 19E,F). . . 11
 Tail not long and thin, but rounded, or conical and more or less pointed
 (Fig. 19, A, B, C, G, H, I) 12
11. Tail tip frequently clavate, (Fig. 19F), one or two ovaries, distance from
 anterior end to center of median oesophageal bulb equal to, or greater
 than distance from center of bulb to base of oesophagus
 Psilenchus
 Tail tip usually pointed, often curved ventrally (Fig. 19E). Distance from
 anterior end to middle of median oesophageal bulb less than distance
 from this point to base of oesophagus. One ovary .. Tylenchus
12. Vulva at 75% or more of body length 13
 Vulva at 60% or less of body length 23
13. Body short and stout, length about 6 to 10 times greatest width 14
 Body slender, length more than 20 times greatest width (Fig. 12K,L). . . 16
14. Body with prominent retrorse (directed posteriorly) annules, usually longer
 than 0.3 mm.; length 10 or more times greatest width. Found in soil . . 15
 Body without prominent annules, length 6 to 8 times greatest width. Length
 about 0.3 mm. Found in roots Cacopaurus
15. Annules with prominent spines or scales Criconema
 Annules without spines or scales (Fig. 12I) Criconemoides
16. Mature female more than 2 mm., often 3-5 mm. long. Found in galls, in
 leaves, or in inflorescence of grains and grasses Anguina
 Body length less than 2 mm. Found in soil or in roots and tubers; sometimes
 in bulbs, leaves and stems 17
17. Stylet longer than 3 times width of lip region (Fig. 18C) 18
 Stylet shorter or about twice width of lip region (Fig. 18A,E) 19

18. Posterior portion of body strongly curved ventrally Paratylenchus
 Body usually nearly straight, usually covered by loose cuticle of fourth
 molt Hemicycliophora
19. Stylet without knobs or with very small knobs; median oesophageal bulb
 occupying nearly full width of body cavity (Fig. 16C, 18E). 20
 Stylet with distinct knobs; median oesophageal bulb occupying less than
 2/3 of width of body cavity (Fig. 16A) 21
20. Tail pointed (Fig. 19G) Aphelenchoides
 Tail rounded (Fig. 19B) Aphelenchus
21. Tail conical, pointed (Fig. 19G) body long and slender; length 40 or more
 times greatest width. Endoparasitic in bulbs, stems, leaves and tubers,
 from mushroom compost, or sometimes in soil Ditylenchus
 Tail tip rounded 22
22. Lip region distinctly set off from body; oesophageal glands forming a
 lobe overlapping intestine; knobs of stylet closely joined. A common
 genus endoparasitic in roots and tubers, also found in soil (Fig. 18A)
 Pratylenchus
 Lip region not distinctly set off from body; oesophageal gland forming a
 distinct basal bulb (Fig. 16A) knobs of stylet separated like an inverted
 Y. A very rare genus from soil Chitinotylenchus
23. At least 2 mm. long, slim; length 45 or more times greatest width. Stylet
 very long, six or more times as long as width of lip region (Fig. 18D).. 24
 Less than 1.5 mm. long, length of stylet not more than 5 times width of
 lip region 25
24. Tail pointed, (Fig. 18C) oesophagus with distinct basal bulb. . Dolichodorus
 Tail rounded, base of oesophagus a lobe overlapping intestine
 Belonolaimus

25. Lip region flattened, stylet length about twice width of lip region Radopholus
 Lip region convex conoid, stylet length 3 or more times width of lip
 region 26
26. Tail 2 or more times as long as anal body diameter, tapering (Fig. 19G).. 27
 Tail shorter than anal body diameter, rounded (Fig. 19A) 28
27. Tail tip rounded Tylenchorhynchus
 Tail tip nearly pointed Tetylenchus
28. Lip region distinctly set off from body, divided into minute plates.
 Body at rest or fixed lying only slightly curved Hoplolaimus
 Annulated lip region continuous with body contour, body usually lies
 in loose spiral when fixed or at rest. (Fig. 12H) 29
29. Opening of dorsal oesophageal gland about one-third stylet length or more
 posterior to stylet knobs Helicotylenchus
 Opening of dorsal oesophageal gland much less than one-third stylet length
 posterior to stylet knobs Rotylenchus

region flattened, stylet length about twice width of lip region

Latipoda

region convex conoid, stylet length 3 or more times width of lip

26

27 or more times as long as anal body diameter, tapering (Fig. 190)

28 than anal body diameter, rounded (Fig. 191)

Latipoda

Tetylema

region distinctly set off from body, divided into minute plates

Hoplostoma only at rest or fixed lying only slightly curved

Annulated lip region continuous with body contour, body usually lies

29 loose spiral when fixed or at rest. (Fig. 192)

Opening of dorsal oesophageal gland about one-third stylet length or more

Latipoda posterior to stylet knobs

Opening of dorsal oesophageal gland much less than one-third stylet length

Botrylenchus posterior to stylet knobs

WHERE NEMATODES
ARE FOUND —




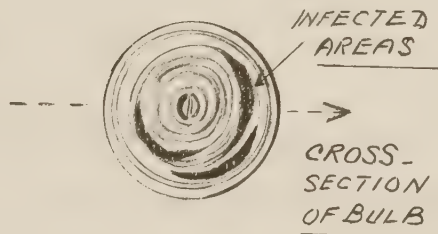



FIGURE - 1.

A.

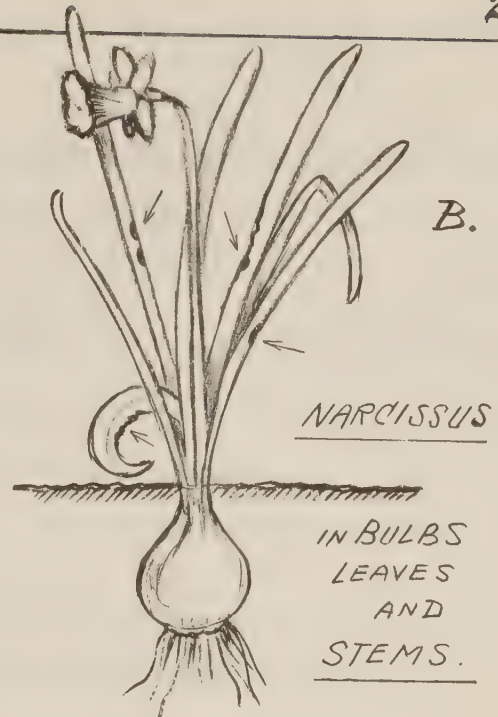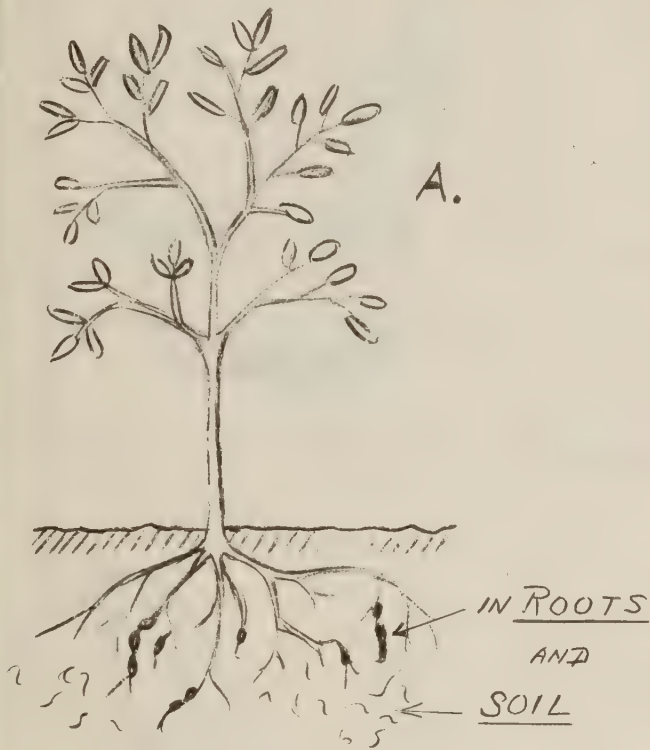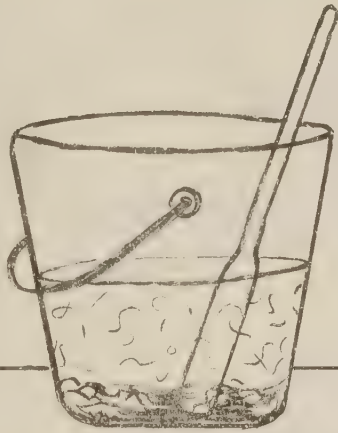

PLACE SOIL SAMPLE IN BUCKET,
ADD WATER AND STIR
WITH STICK.
SAND SETTLES TO BOTTOM
NEMATODES REMAIN
SUSPENDED

B.

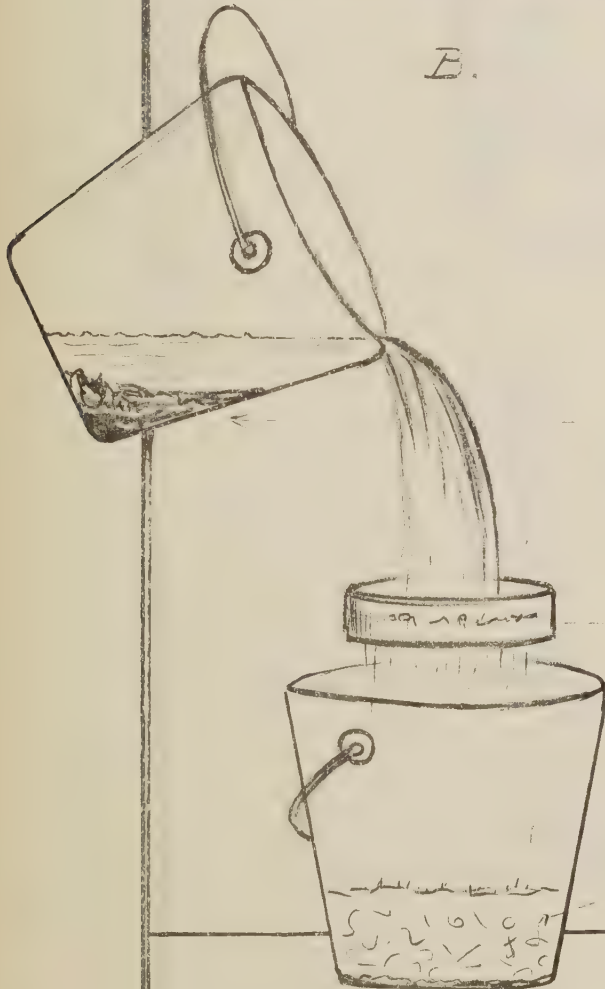


POUR THROUGH 20-MESH SIEVE
INTO SECOND BUCKET.

SAND AND GRAVEL REMAIN
IN FIRST BUCKET

SIEVE REMOVES LARGE
PIECES OF DEBRIS

NEMATODES PASS
THROUGH SIEVE

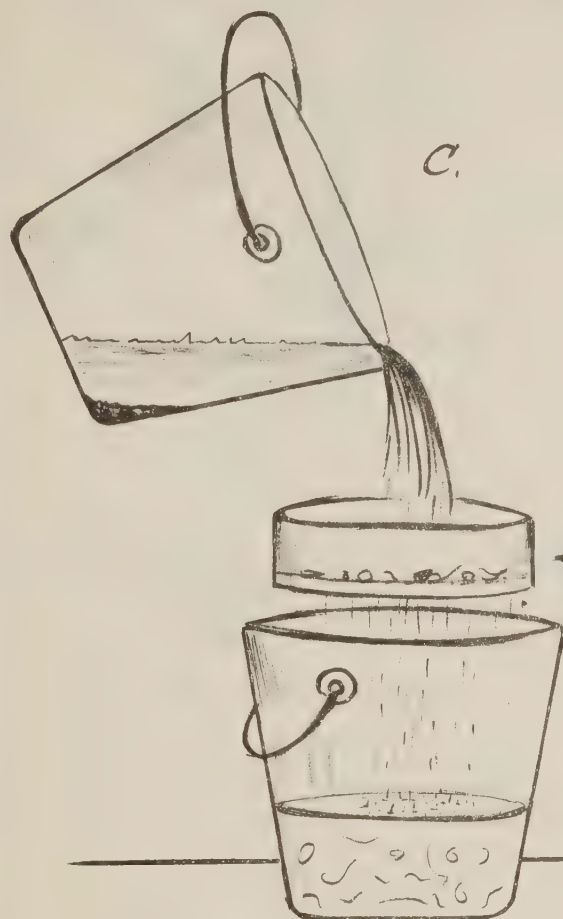


POUR THROUGH 60-MESH SIEVE

MORE SAND SETTLES OUT

← CYSTS AND LARGE
NEMATODES
REMAIN ON SIEVE.

THE SMALLER
NEMATODES IN WATER.

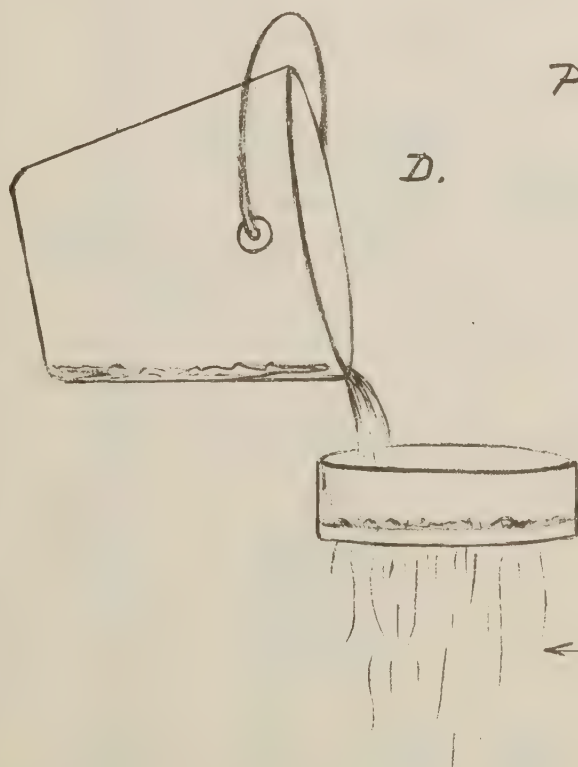

POUR THROUGH 200-MESH

SMALL NEMATODES
REMAIN ON SIEVE
WITH PLANT DEBRIS

← WATER AND FINE SOIL
PARTICLES PASS THROUGH.

FIGURE - 3.

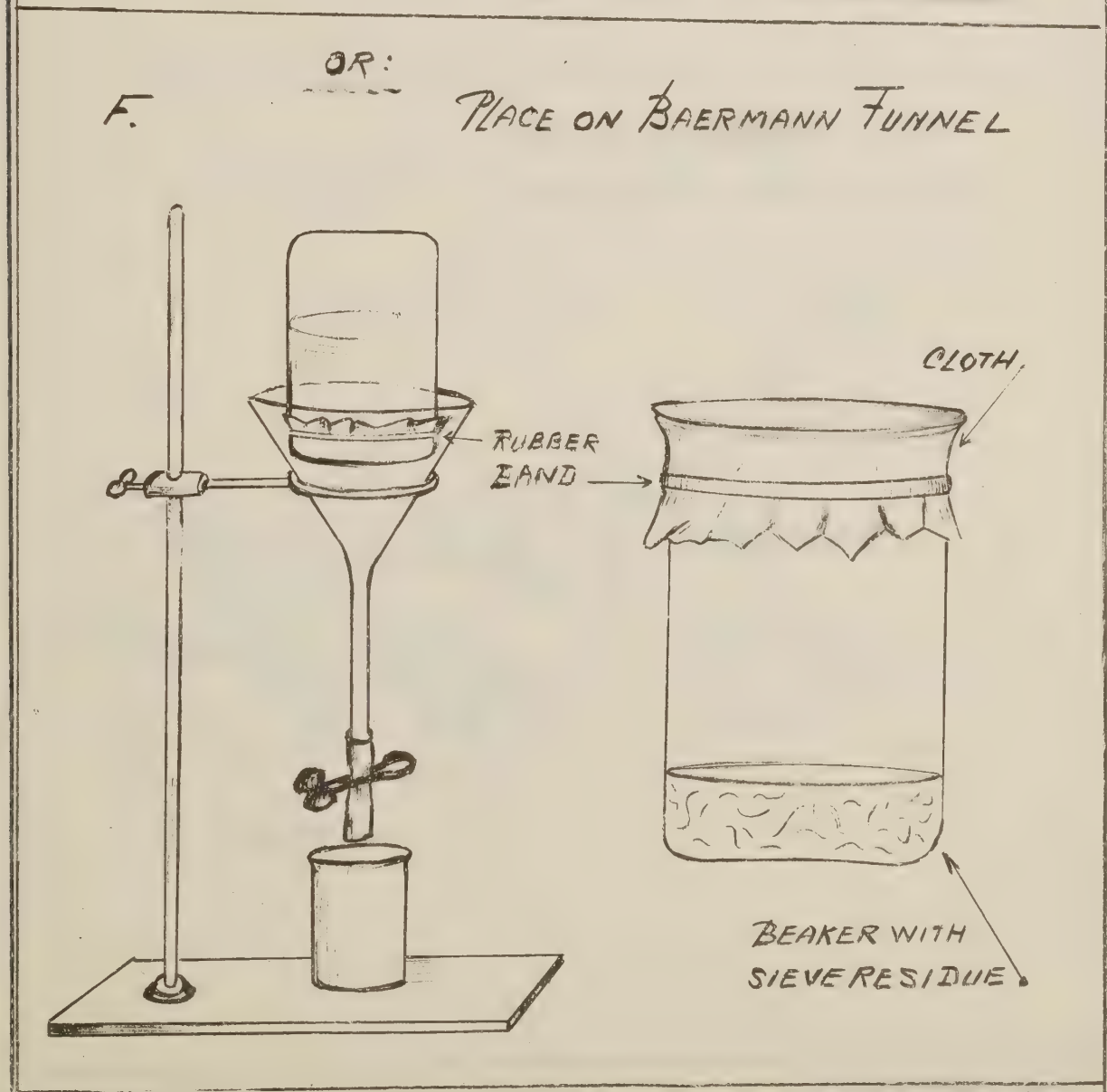
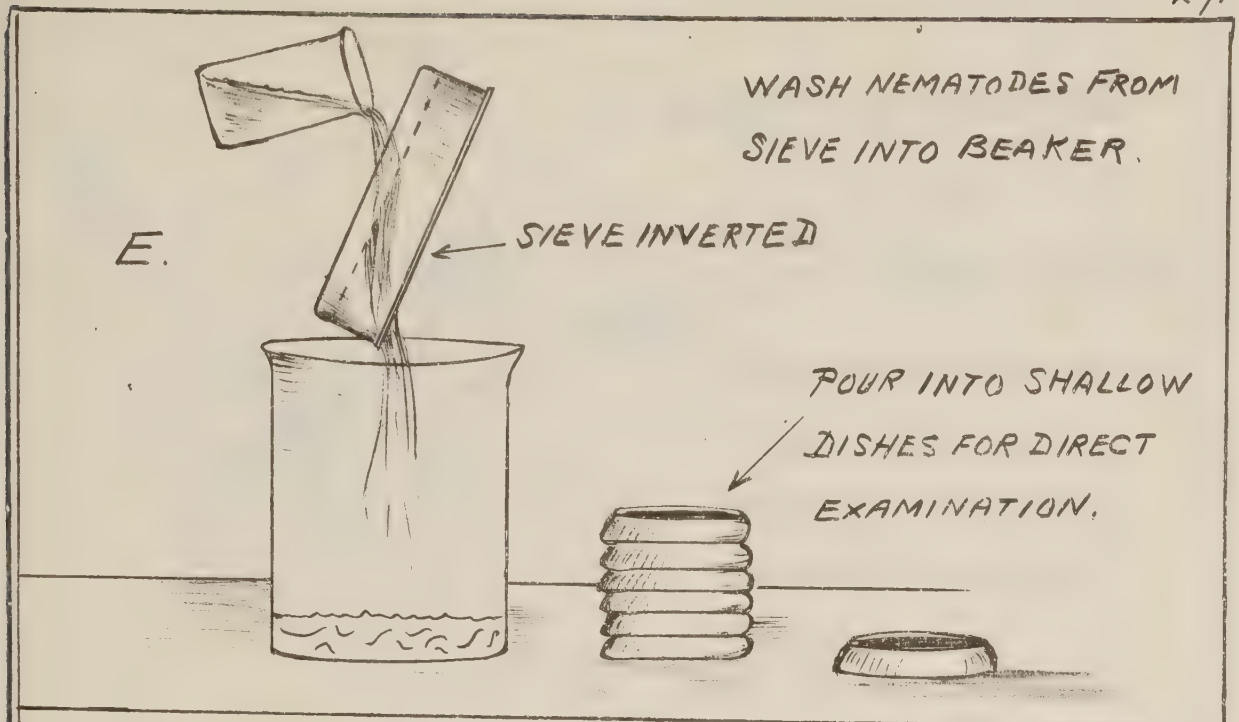


FIGURE - 4.




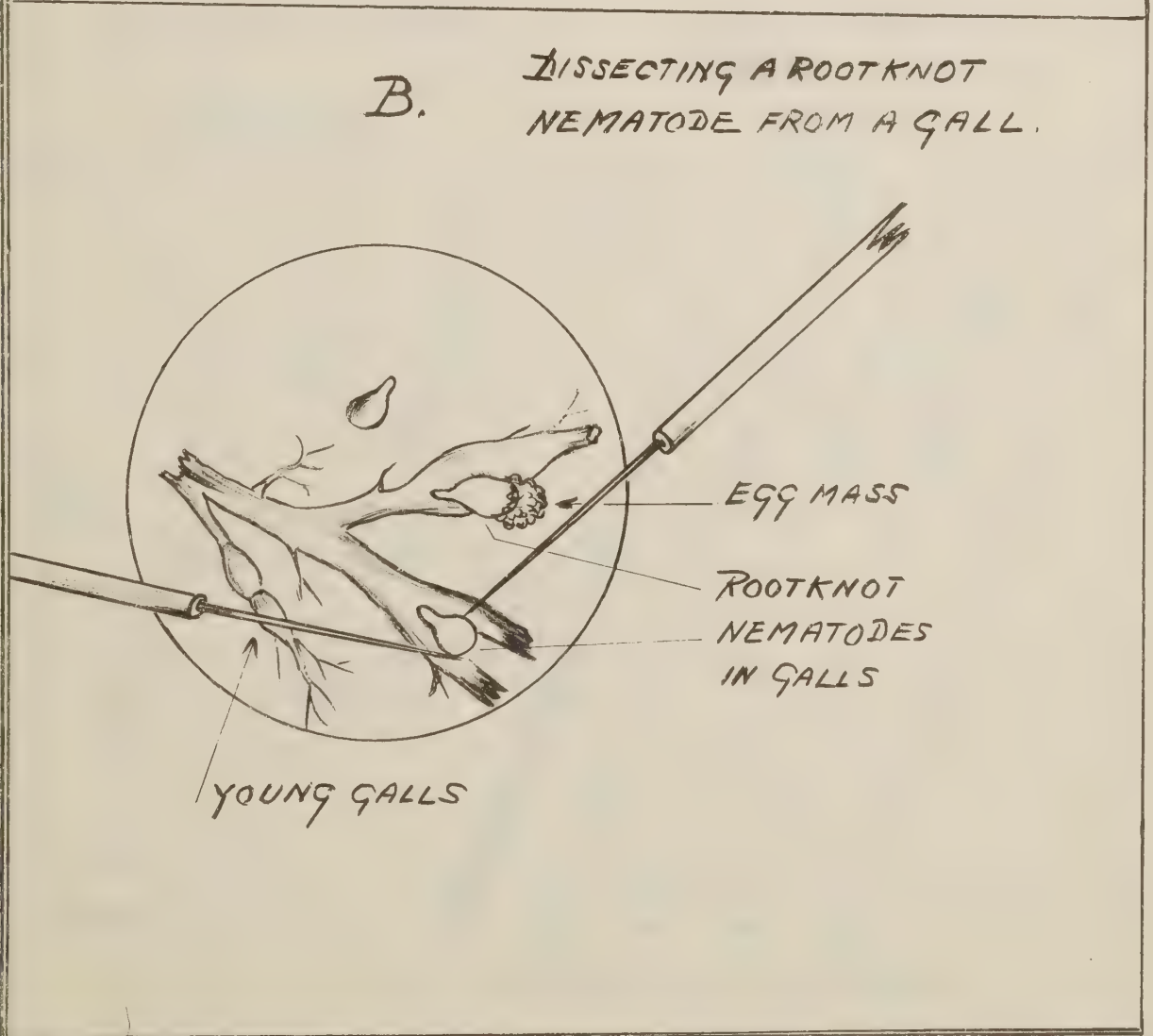


FIGURE - 5



A.

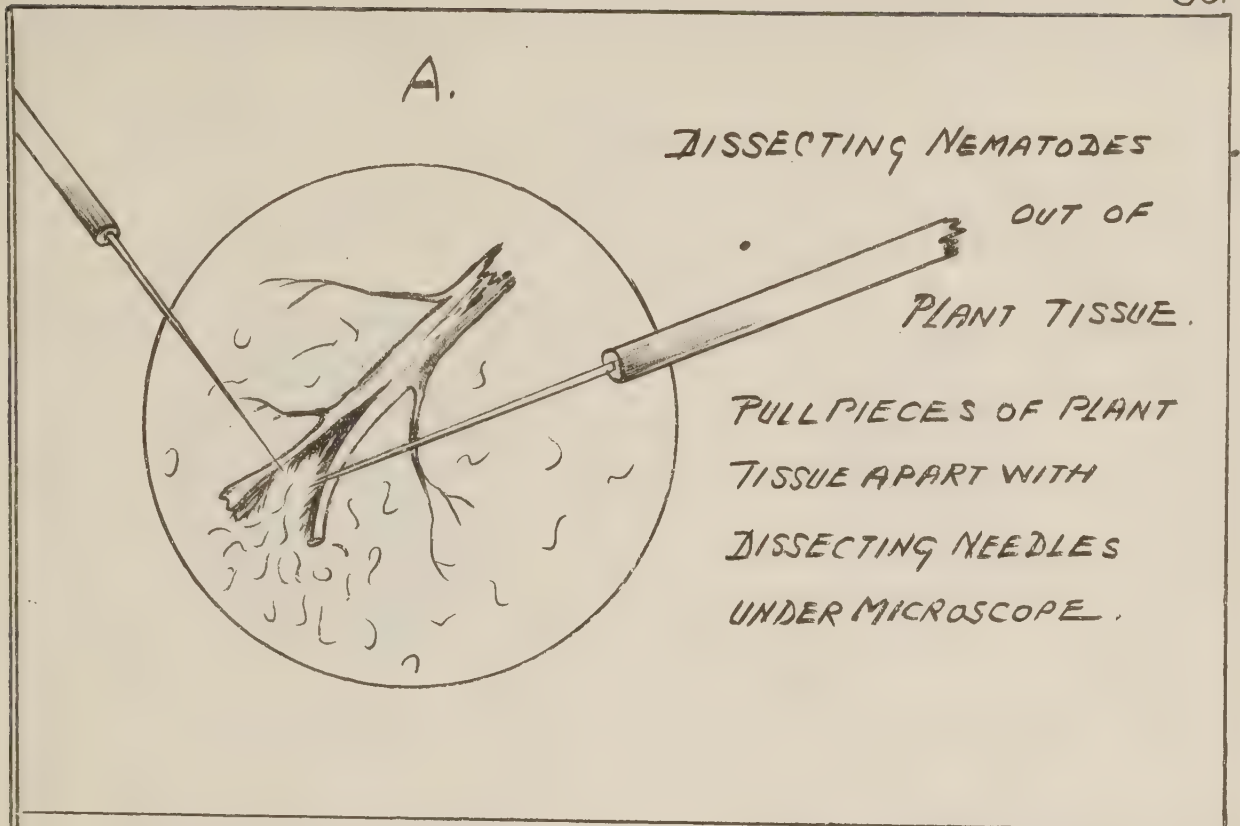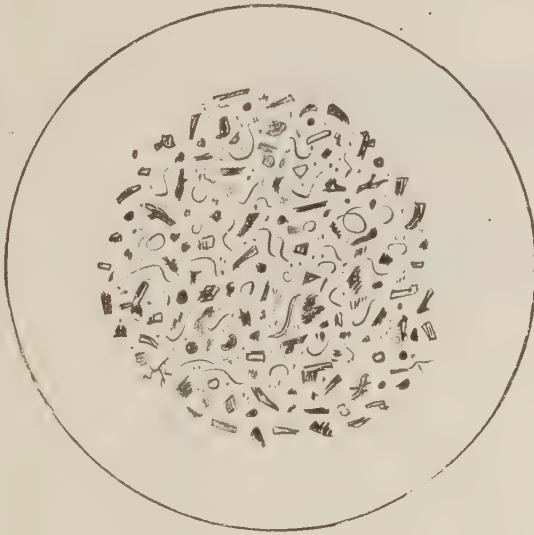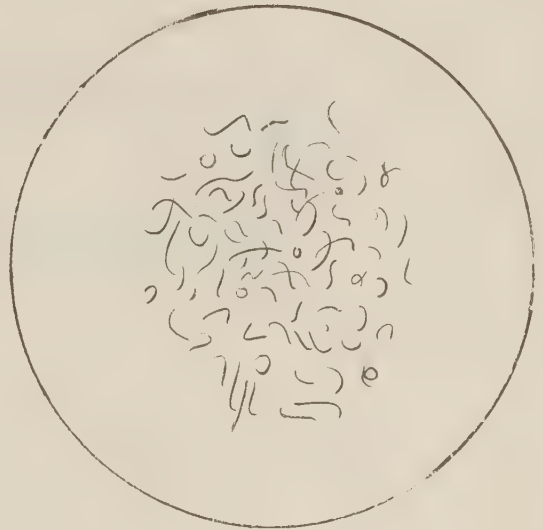


NEMATODES IN DEBRIS.

OBSCUREDNEMATODES FROM
BAERMANN FUNNELCLEAN

B.

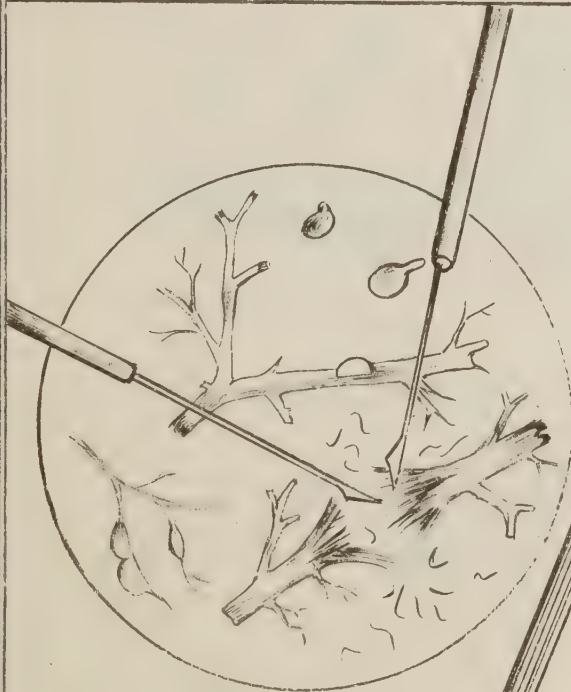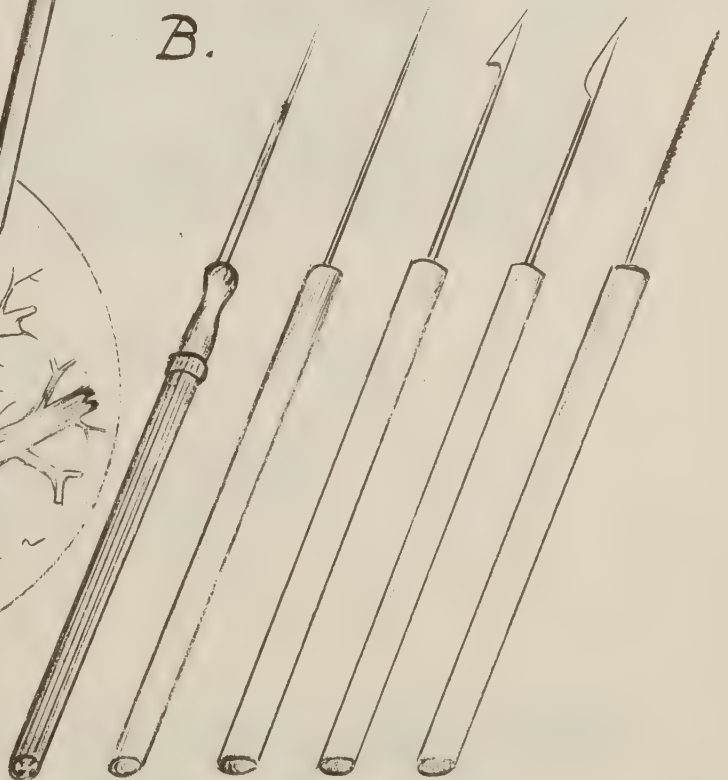DISSECTING
PLANT ROOTS

DISSECTING NEEDLES

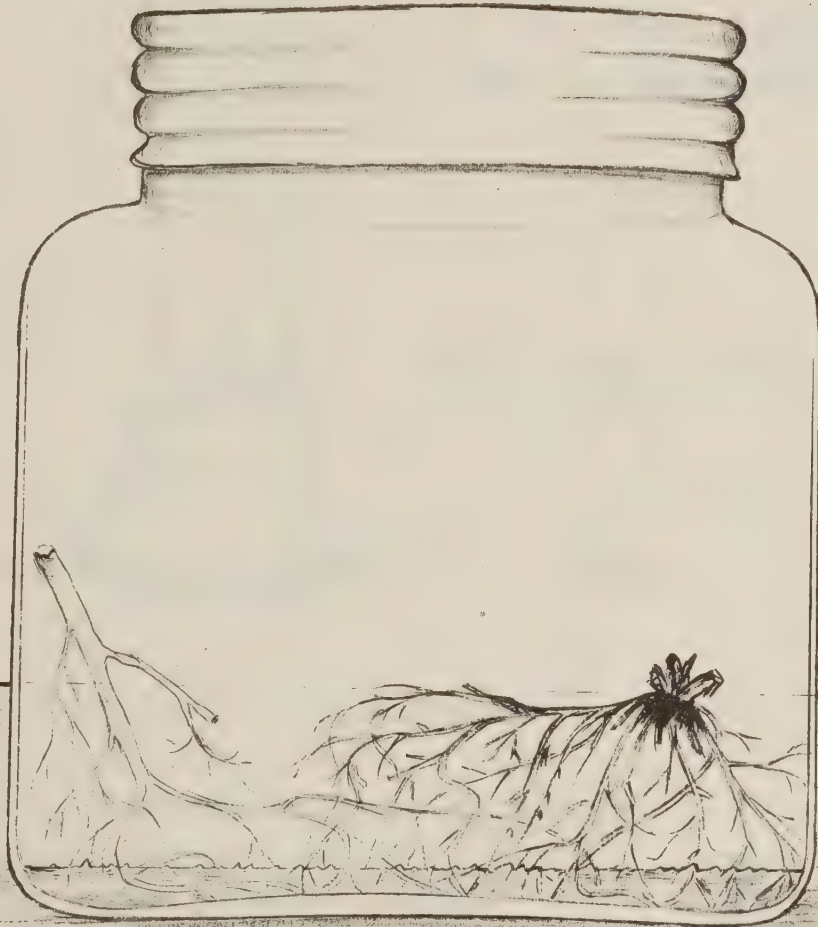

INCUBATION TECHNIQUE
FOR
EXTRACTING NEMATODES
FROM ROOTS


WET ROOTS IN FRUIT JAR; ADD A LITTLE WATER;
POUR OFF WATER AND EXAMINE FOR NEMATODES.

WARING BLENDOR

1. PLACE PLANT MATERIAL
IN BLENDOR WITH WATER

2. OPERATE BLENDOR
FOR 30 SECONDS

A.

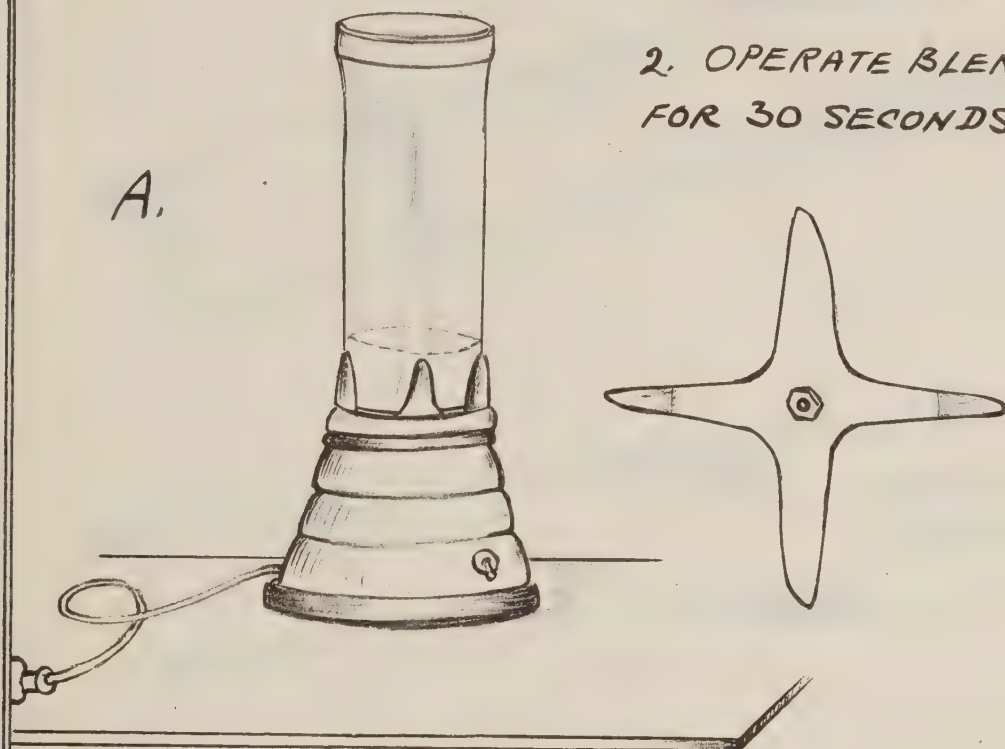

B.

3. PLACE MATERIAL FROM BLENDOR
IN 60-MESH SIEVE
OVER 200-MESH SIEVE

4. WASH WITH TAP WATER

5. WASH NEMATODES
FROM 200-MESH
SIEVE AS SHOWN
IN FIGURE 4.E.



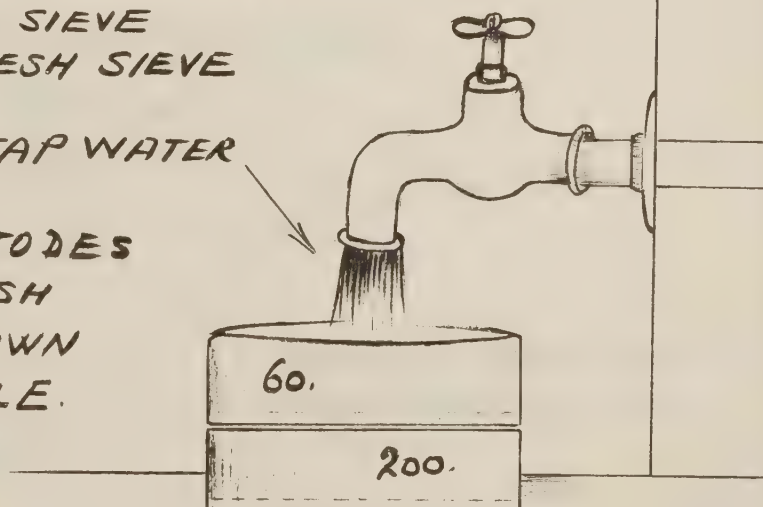



FIGURE - 8.

HARD BAMBOO

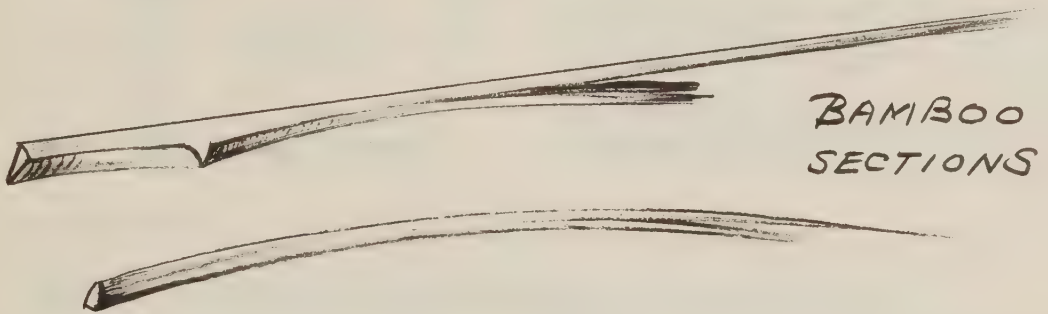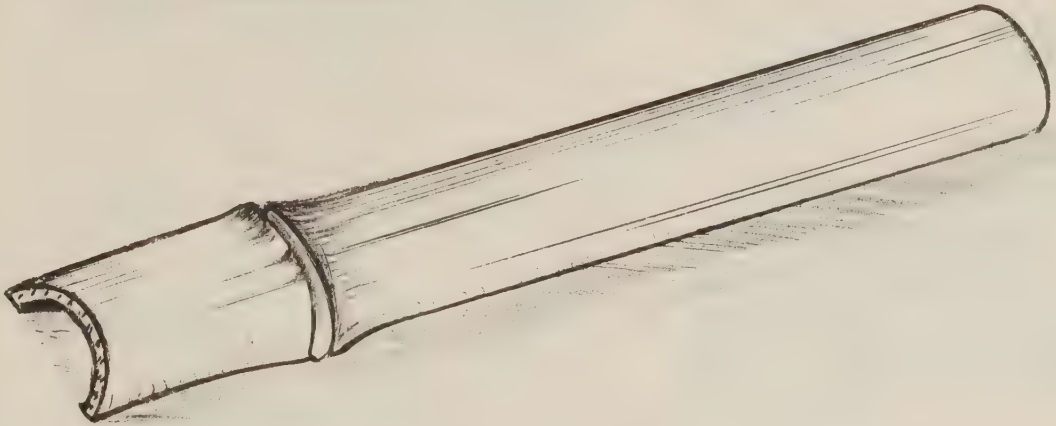

BAMBOO
SECTIONS

BAMBOO SPLINTER
TRIMMED AND
POINTED

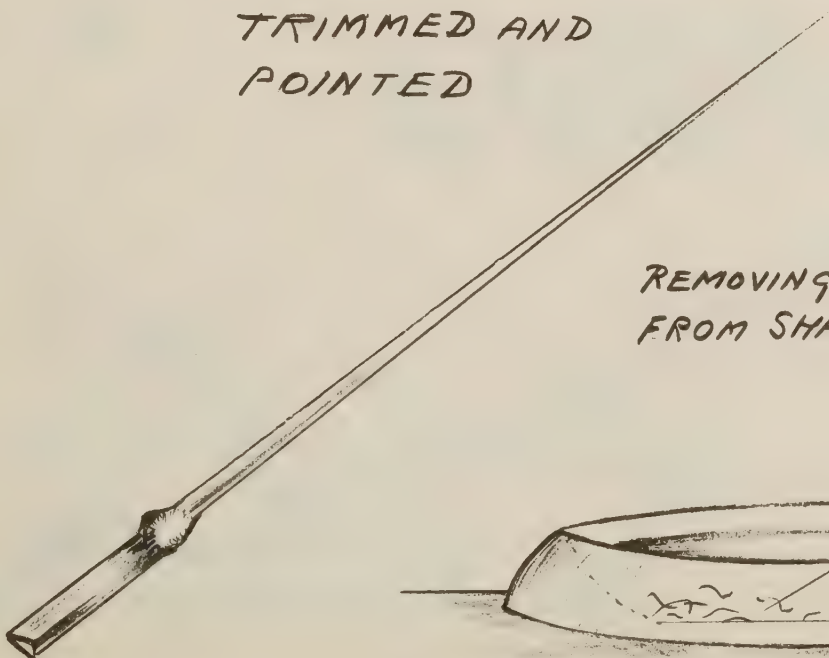

REMOVING NEMATODES
FROM SHALLOW DISH.


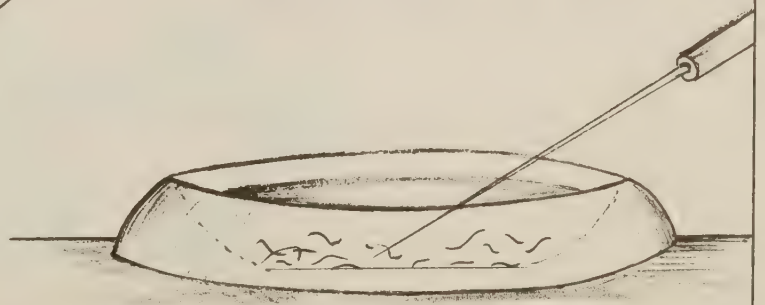



FIGURE - 9.

MOUNTING NEMATODES ON MICROSCOPE SLIDES

A.

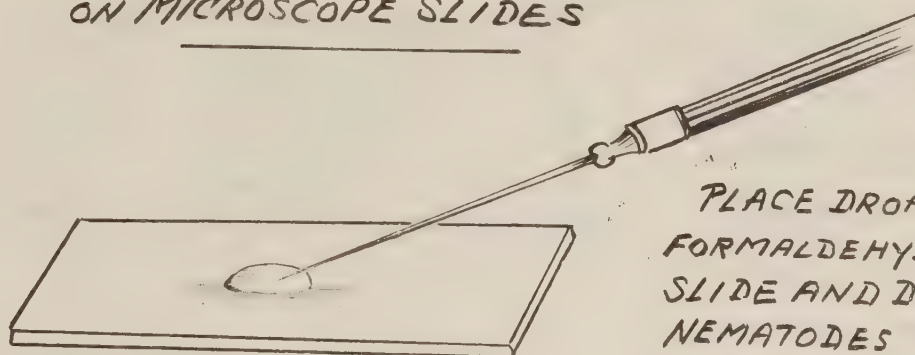

PLACE DROP OF 5%
FORMALDEHYDE ON
SLIDE AND DEPOSIT
NEMATODES IN THE
DROP, USING BAMBOO
SPLINTER.

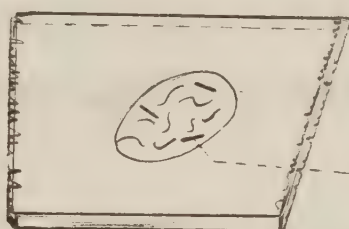

YOU CAN MOUNT AS MANY AS
20 NEMATODES ON ONE SLIDE.

PIECES OF GLASS ROD TO SUPPORT
COVER GLASS (IF DESIRED)

LOOK AT THE DROP UNDER THE DISSECTING
MICROSCOPE AND BE SURE THAT THE NEMATODES
LIE IN THE BOTTOM OF THE DROP AND
IN CONTACT WITH THE SLIDE.

B.

USING FORCEPS,
PLACE A COVER GLASS
OVER THE DROP.



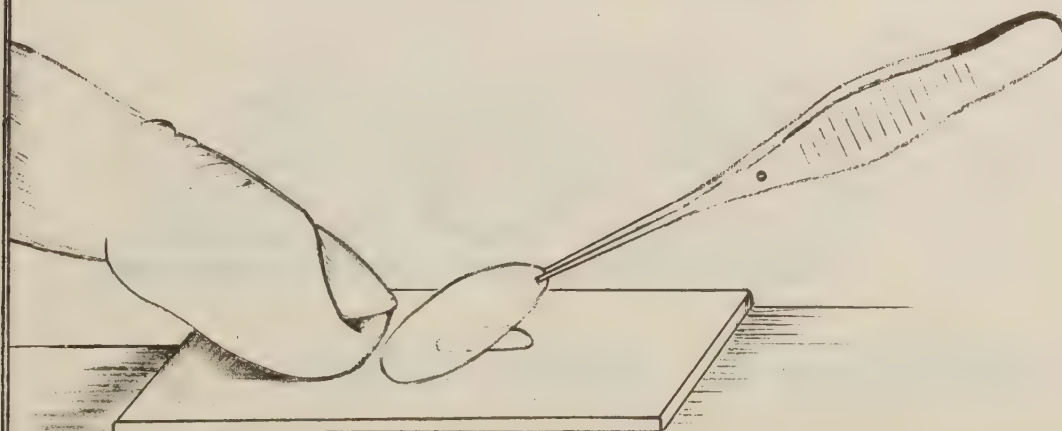


FIGURE-10.

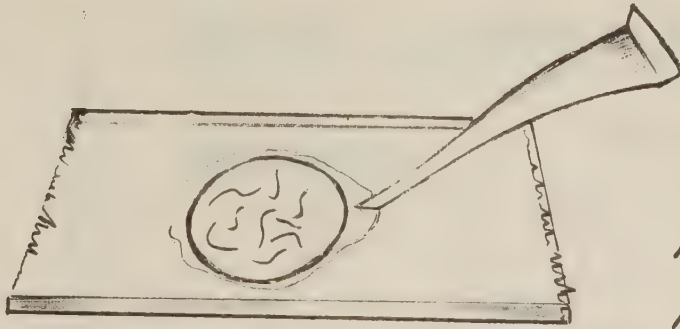


USE A SMALL
PIECE OF FILTER
PAPER, BLOTTING
PAPER OR PAPER
TOWEL TO ABSORB
EXCESS WATER

DO THIS UNDER THE
DISSECTING MICROSCOPE TO BE
SURE THAT NO NEMATODES ARE LOST. A.

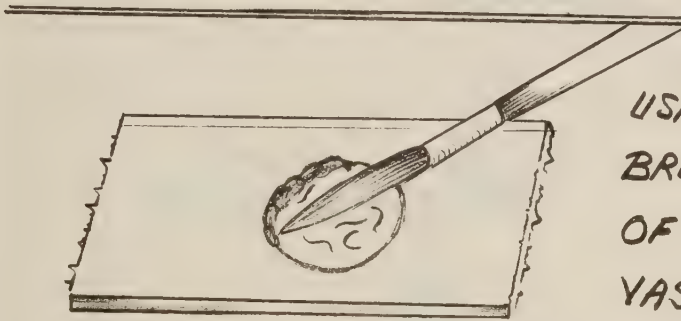

USING A SMALL B.
BRUSH APPLY A RING
OF HOT PARAFFIN-
VASELINE MIXTURE,
ZUT, OR OTHER
SEALING MATERIAL.

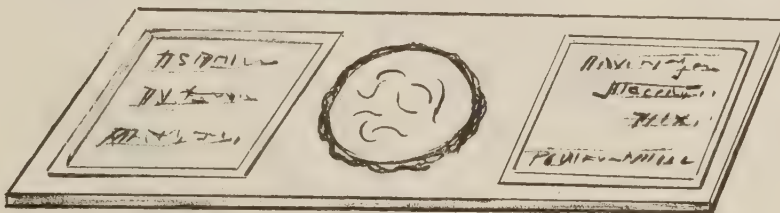

C.
LABEL SLIDE, USING FINE PEN.

SHAPES OF VARIOUS NEMATODES

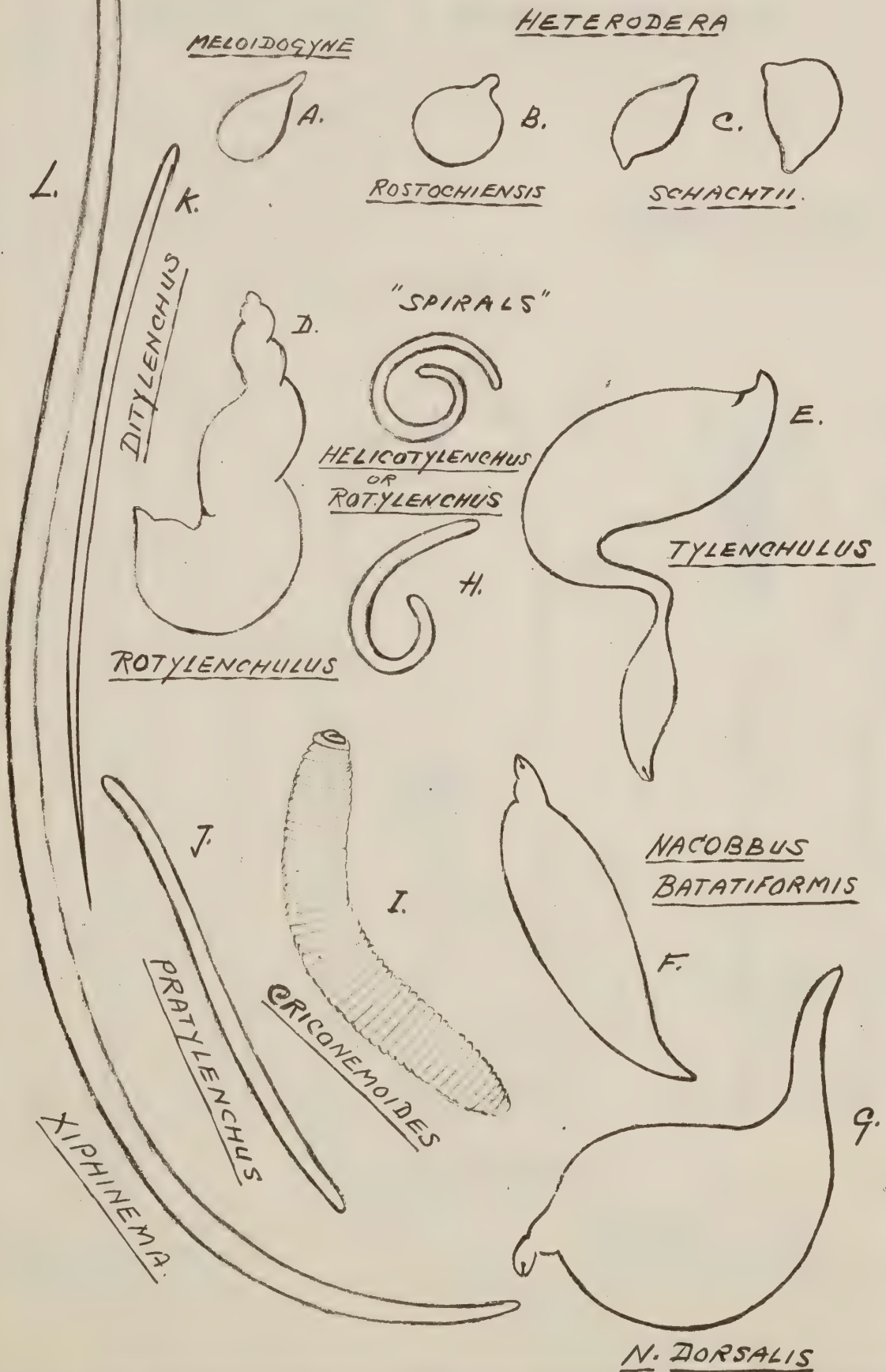



FIGURE 12.

VARIATIONS' in BUCCAL CAPSULES.


TUBULAR

ODONTOSTYLET.


STOMATOSTYLET.

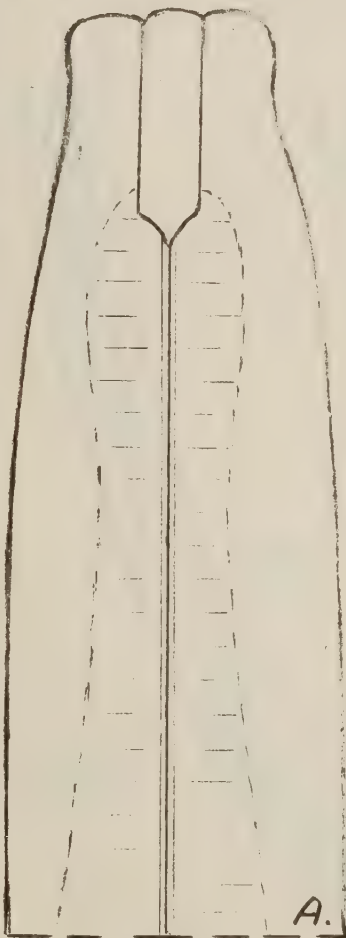

A.

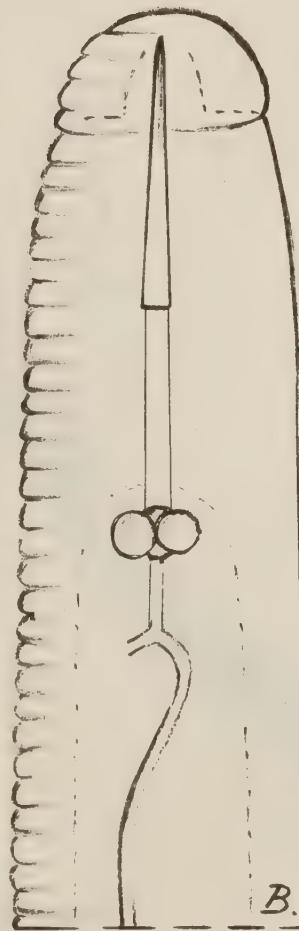

B.

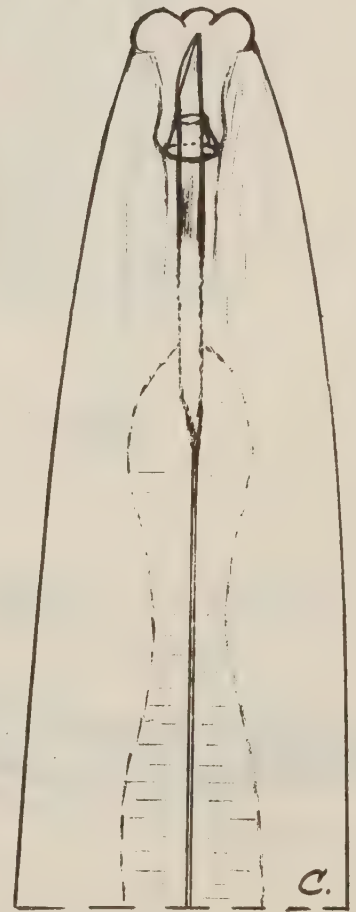

C.



RHABDITIS

ROTYLENCHUS

DORYLAIMUS


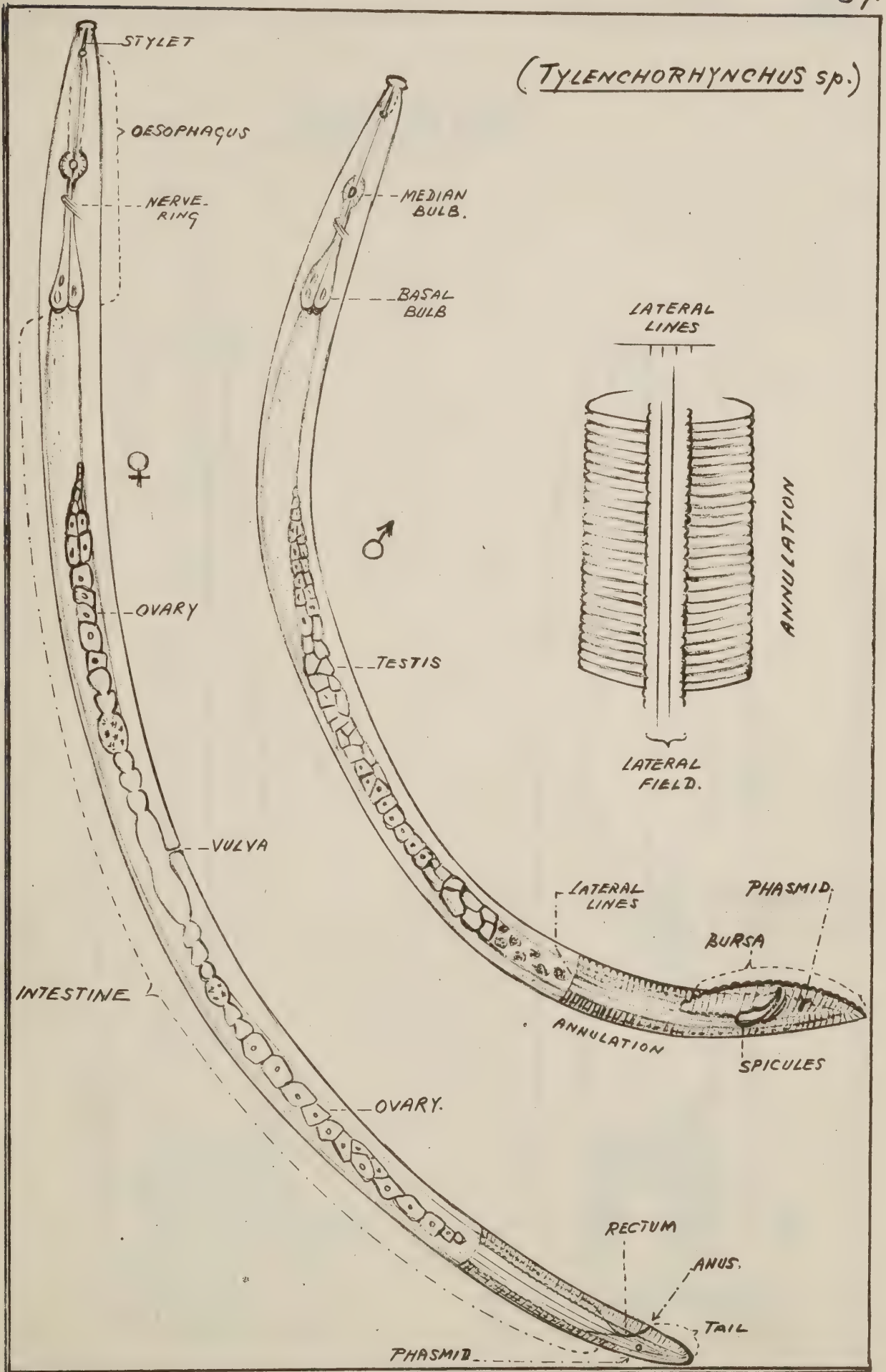



FIGURE 14.

PHASMIDIA.

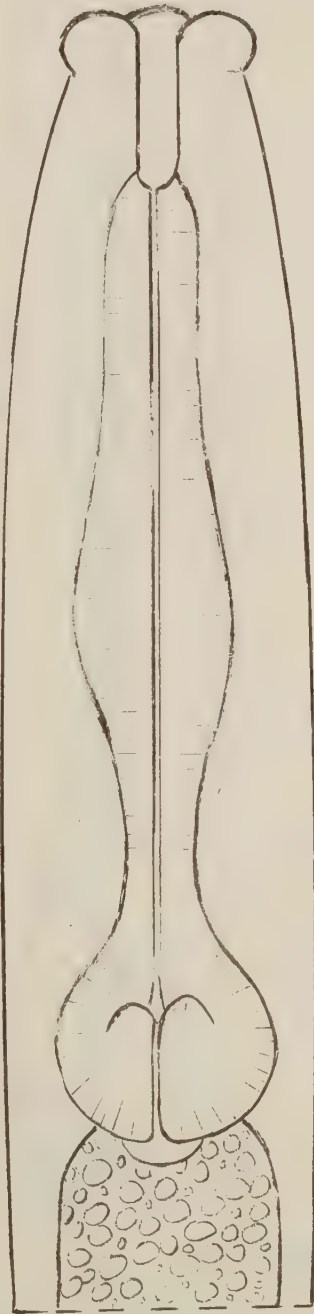

A.

RHABDITOID

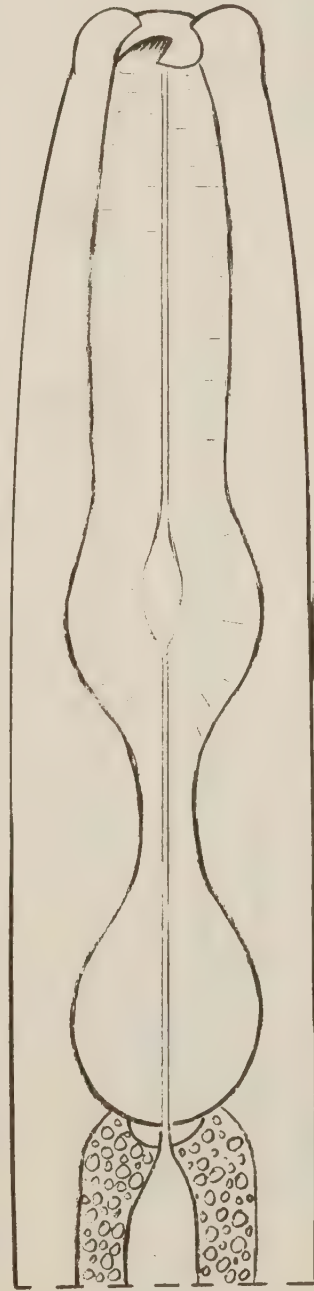

B.

DIPLOGASTEROID.


PHASMIDIA



A.



B.

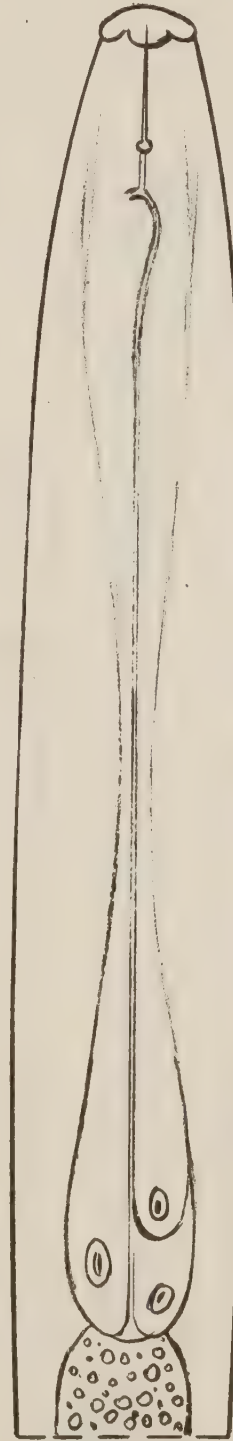


C.

APHELENCHOID


TYLENCHOID

NEOTYLENCHOID


APHASMIDIA.

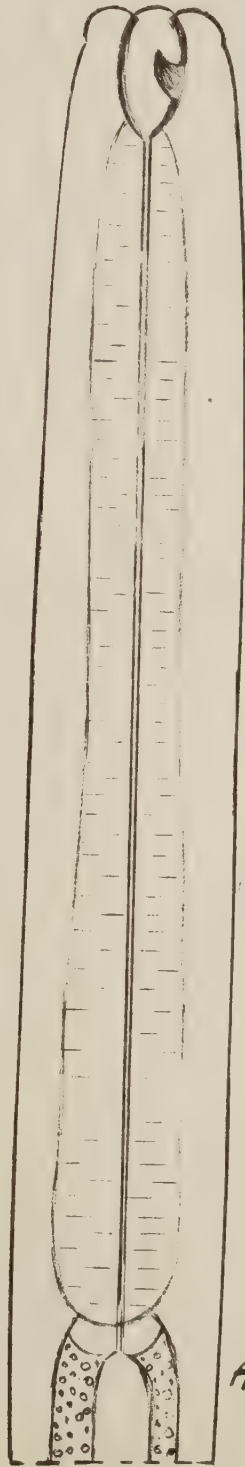

A.

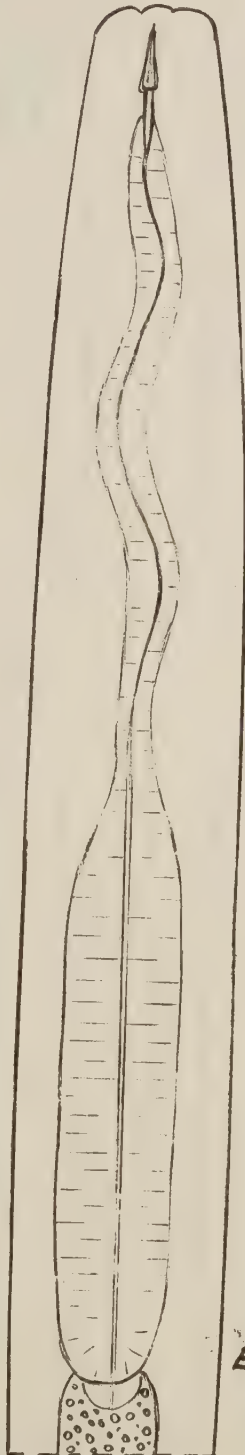CYLINDROID.


B.

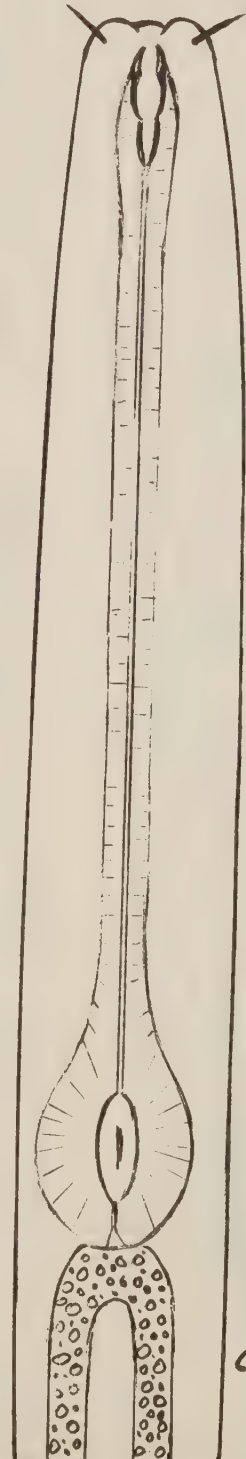DORYLAIMOID.


C.

PLECTOID.


STYLET SHAPES

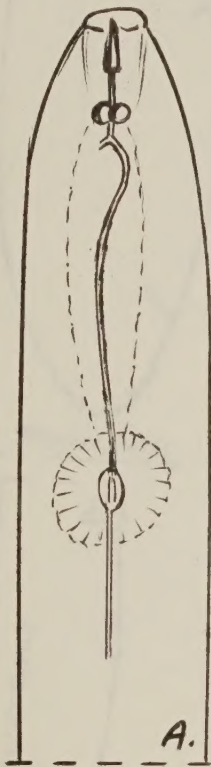


PRATYLENCHUS

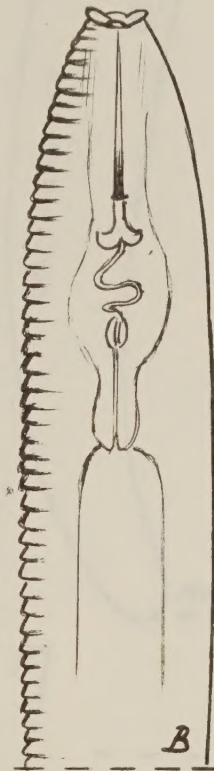


CRICONEMOIDES

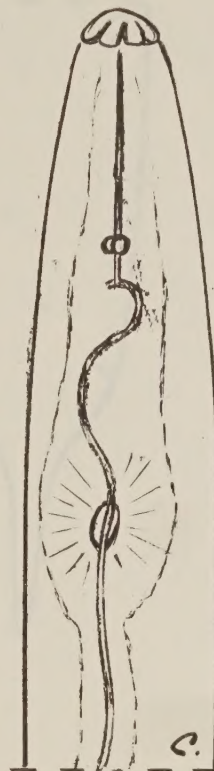


PARATYLENCHUS

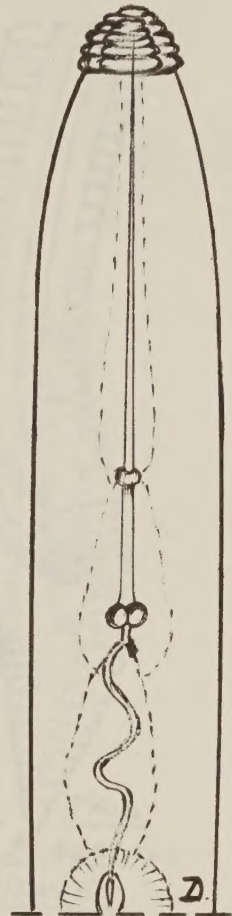


BELONDOLAIMUS

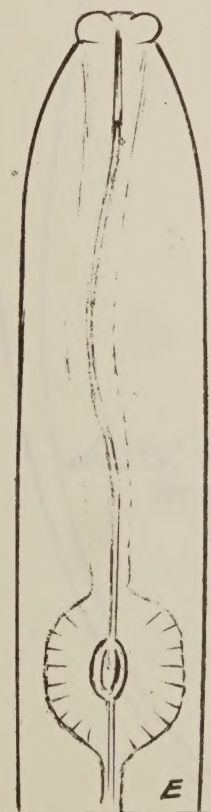


APHELENCHOIDES

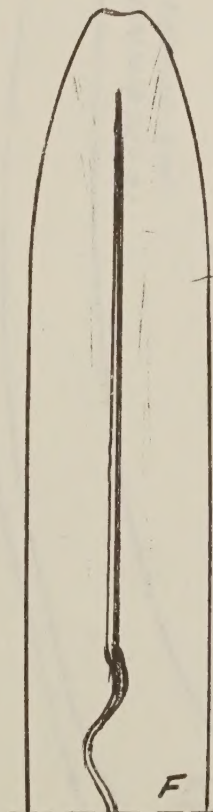

LONGIDORUS

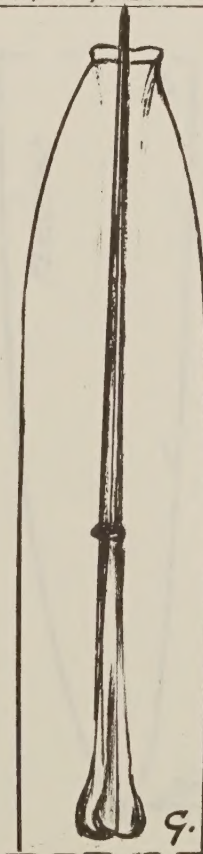


XIPHINEMA


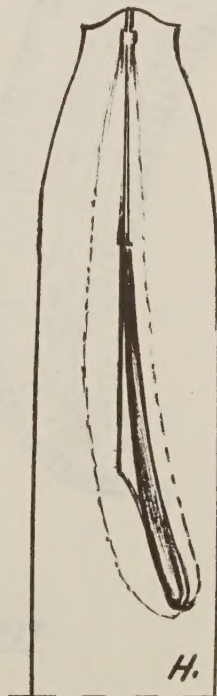


TRICHODORUS


TAIL SHAPES

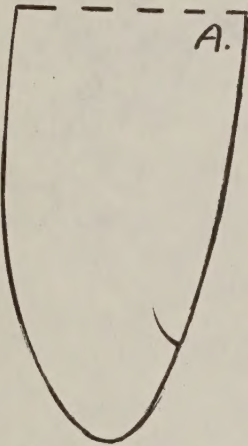





ROUNDED

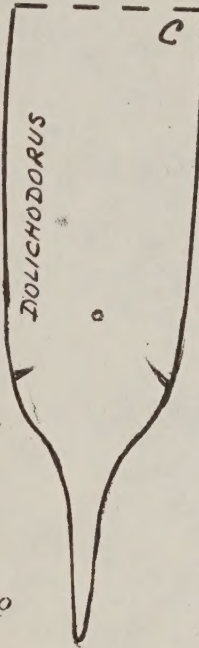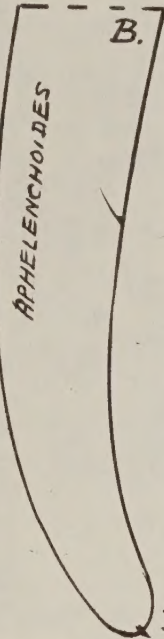

MUCRO

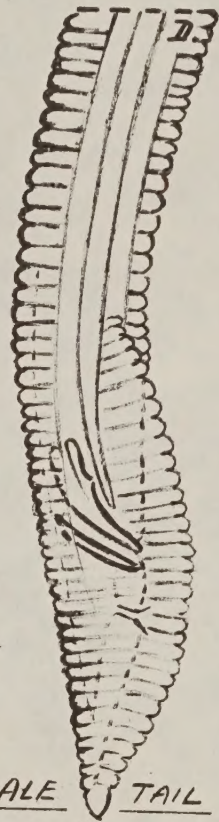

MALE TAIL
WITH BURSA AND
SPICULES

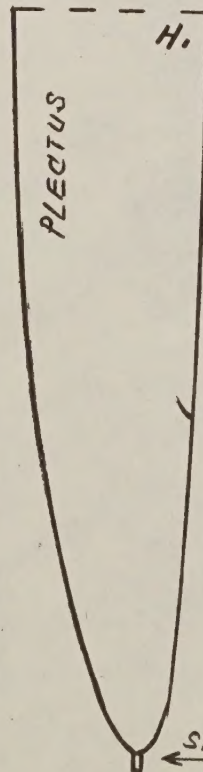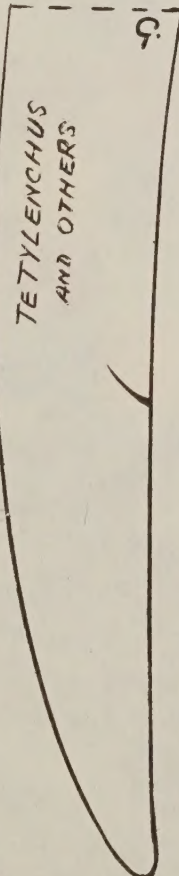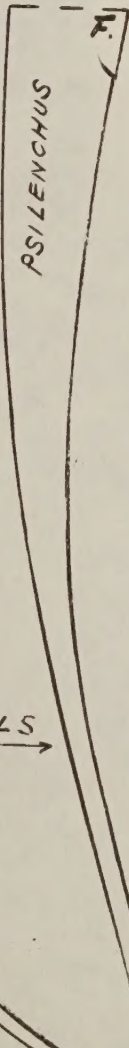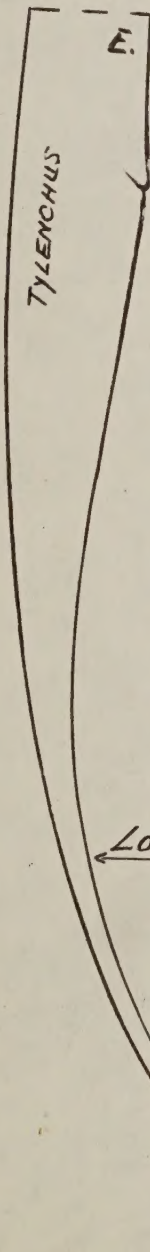

SPINNERET

POINTED

LONG
TAILS


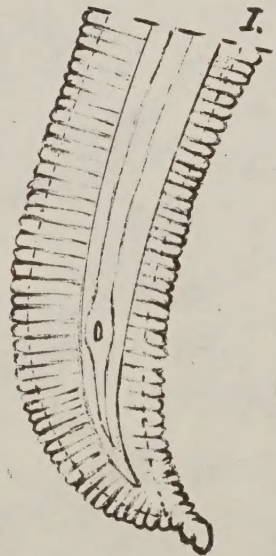


FIGURE 19.







1022507867


TAIL SHAPES

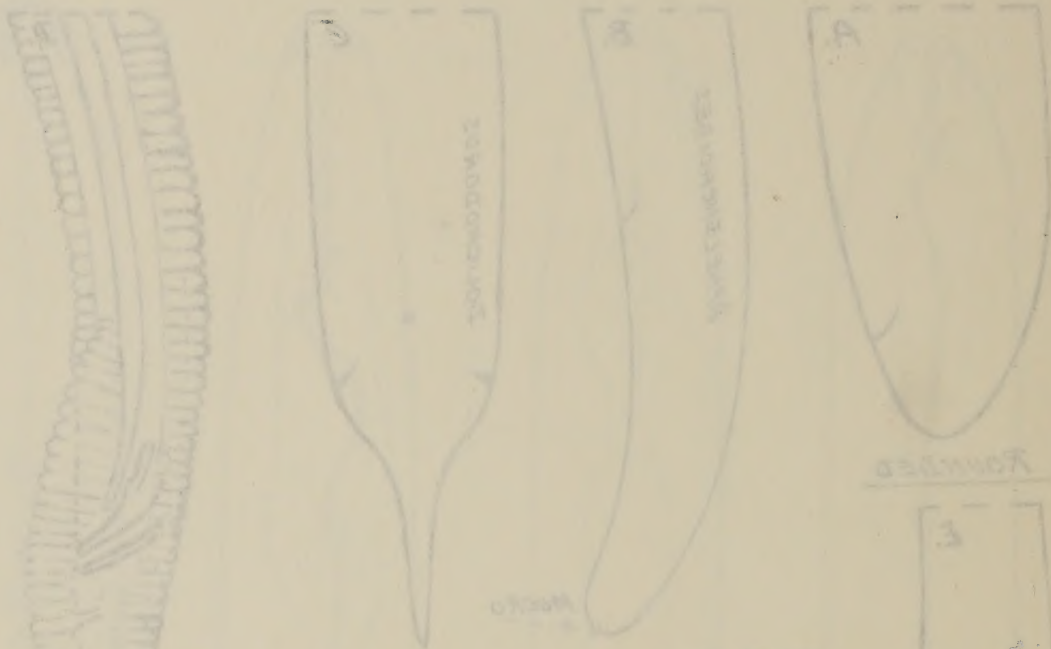

ROUNDED

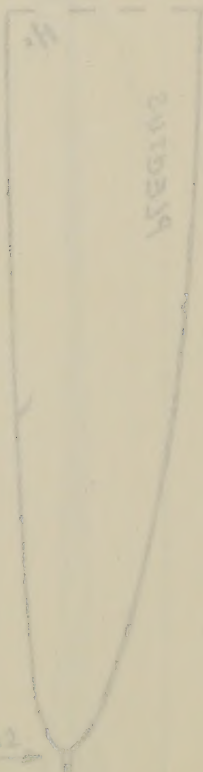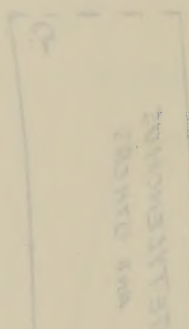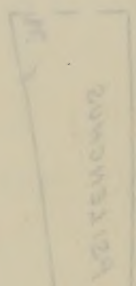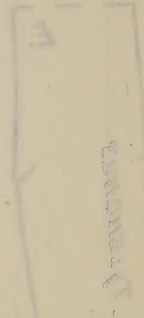

MALE
WITH BARS AND
SPINES

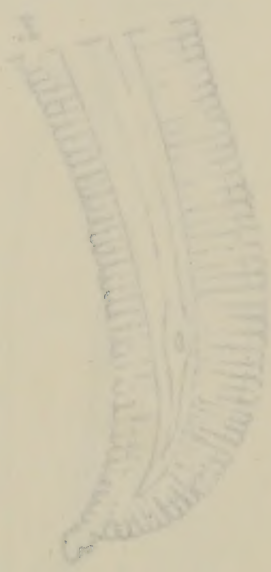

SPINNY

POINTED

LONG
TAILS